


FFAP Test Mix**Column:** Zebron[™] ZB-FFAP, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G009-11**Oven Profile:** 140°C Isothermal**Carrier Gas:** Hydrogen, 1.93 mL/min**Injection:** Split 50:1 1 µL @ 250°C**Detection:** Flame Ionization (FID) (255°C)**Products used in this application:**

FFAP Test Mix

ANALYTES:

- 1 Methane
- 2 Tridecane
- 3 Tetradecane
- 4 1-Octanol
- 5 Methyl decanoate
- 6 Isovaleric acid
- 7 Methyl undecanoate
- 8 Naphthalene
- 9 Methyl dodecanoate
- 10 1-Hexanoic acid
- 11 2,6-Dimethylphenol

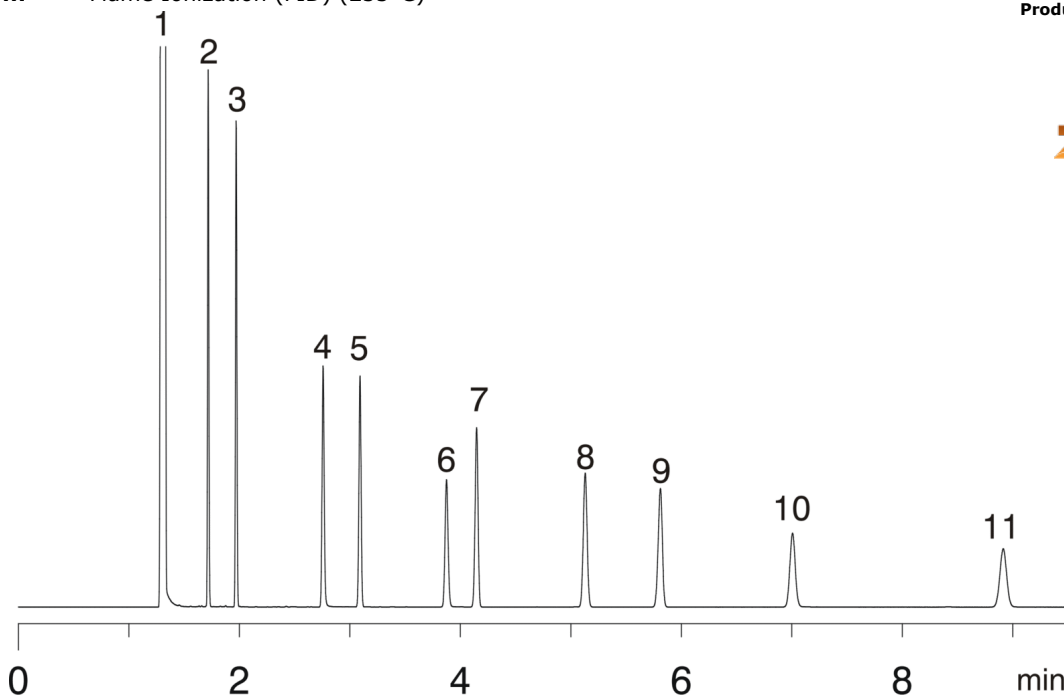