

# HPLC Application

ID No.: 10851

## Up to 60% Reduction in Analysis Time

**Column:** Luna® 3 µm C18(2) 100 Å, LC Column 20 x 4 mm, Ea

**Dimensions:** 20 x 4 mm ID

**Order No:** 00M-4251-D0

**Elution Type:** Gradient

**Eluent A:** Water with 0.1% formic acid

**Eluent B:** Acetonitrile with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	2.4	5	95

**Flow Rate:** 1.6 mL/min

**Col. Temp.:** ambient

**Detection:** Refractive Index @ 0.0000000000 (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

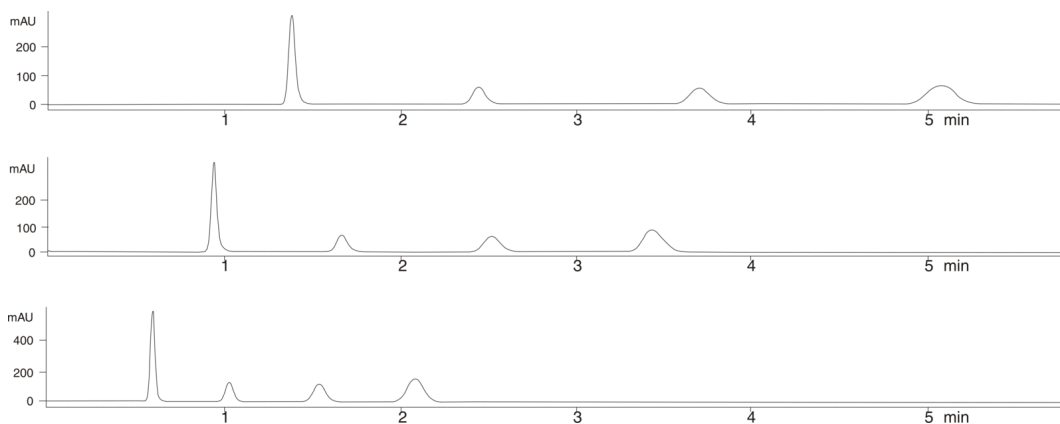

### ANALYTES:

- 1 Acetaminophen
- 2 Propranolol
- 3 Imipramine
- 4 Naproxen
- 5 Valerophenone

