HPLC Application

ID No.: 10872



Increase Sample Through-Put (fig. 2)

Luna® 3 µm C18(2) 100 Å, LC Column 20 x 4 mm, Ea

20 x 4 mm ID **Dimensions:** Order No: 00M-4251-D0 Elution Type: Isocratic

Eluent A: Water with 0.1% ammonium hydroxide, pH 10 Acetonitrile with 0.1% ammonium hydroxide, pH 10 Eluent B: Gradient Step No. Time (min) Pct A Pct B **Profile:** 80 20

2 mL/min Flow Rate: Col. Temp.: 25 °C

UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient) **Detection: Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:



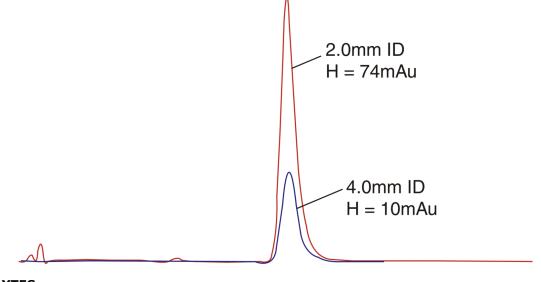

ANALYTES:

1 Pindolol

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at support@phenomenex.com

