

Basic Drugs using Ammonium Hydroxide

Column: Luna® 3 µm C18(2) 100 Å, LC Column 30 x 4.6 mm, Ea

Dimensions: 30 x 4.6 mm ID

Order No: 00A-4251-E0

Elution Type: Gradient

Eluent A: Water with 0.1% ammonium hydroxide, pH 10

Eluent B: Acetonitrile with 0.1% ammonium hydroxide, pH 10

Gradient	Step No.	Time (min)	Pct A	Pct B
----------	----------	------------	-------	-------

Profile:	1	0	85	15
----------	---	---	----	----

2 3 25 75

Flow Rate: 2 mL/min

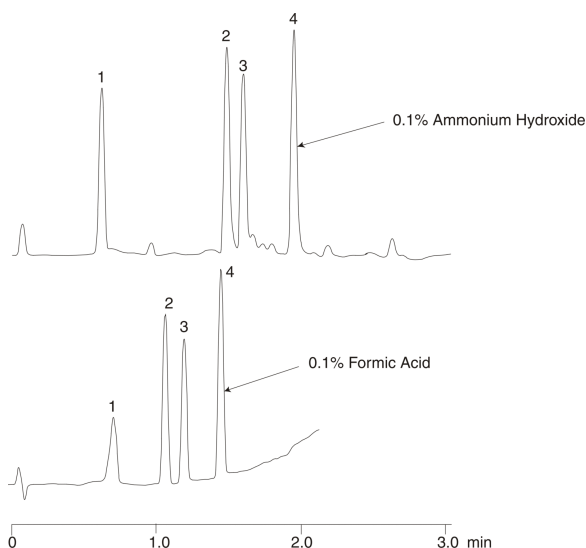
Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 205 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Metaproterenol
- 2 Pindolol
- 3 Metoprolol
- 4 Propranolol



Products used in this application:

