

## Steroids

**Column:** Luna® 3 µm C18(2) 100 Å, LC Column 75 x 4.6 mm, Ea

**Dimensions:** 75 x 4.6 mm ID

**Order No:** 00C-4251-E0

**Elution Type:** Isocratic

**Eluent A:** 0.1 H<sub>3</sub>PO<sub>4</sub>

**Eluent B:** Acetonitrile

**Eluent C:** Methanol

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B | Pct C |
|-------------------|----------|------------|-------|-------|-------|
| 1                 | 1        | 0          | 54    | 35    | 11    |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

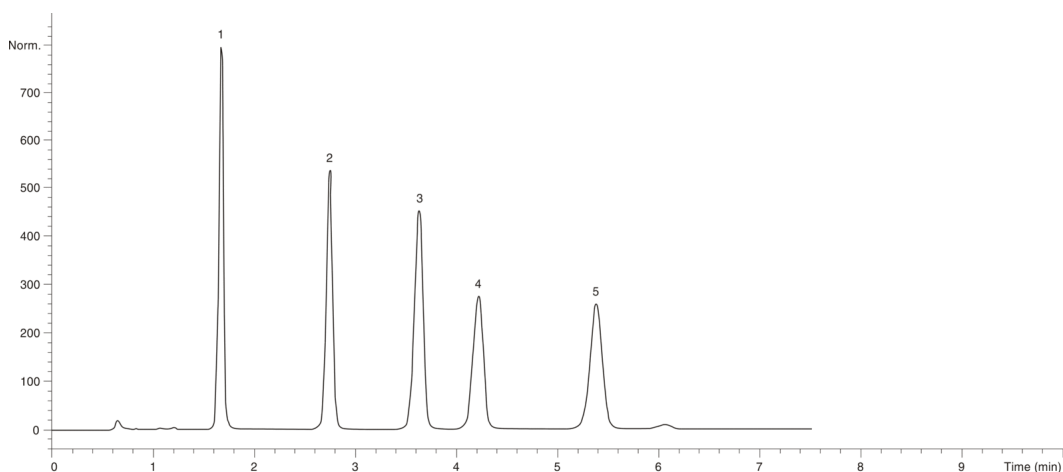

## ANALYTES:

- 1 Hydrocortisone
- 2 Corticosterone
- 3 11- $\alpha$ -Hydroxyprogesterone
- 4 Cortisone acetate
- 5 11-Ketoprogesterone

