


Solvents and Halogenated Hydrocarbons**Column:** Zebron[™] ZB-35, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G003-11**Oven Profile:** 32-125°C @ 12°C (hold 1min)**Carrier Gas:** Constant Flow Helium, 1.12 mL/min**Injection:** Split 50:1 1 μ L @ 250°C**Detection:** Mass Spectrometer (MS) (250°C)**Products used in this application:**

Solvents and Halogenated Hydrocarbons

ANALYTES:

- 1 trans-1,2-Dichloroethylene
- 2 1,1,1-Trichloroethane
- 3 Benzene
- 4 1,2-Dichloroethane
- 5 Bromodichloromethane
- 6 cis-1,3-Dichloropropene
- 7 Toluene
- 8 trans-1,3-Dichloropropene
- 9 Ethylbenzene
- 10 Bromoform
- 11 1,1,2,2-Tetrachloroethane
- 12 1,3-Dichlorobenzene
- 13 1,4-Dichlorobenzene

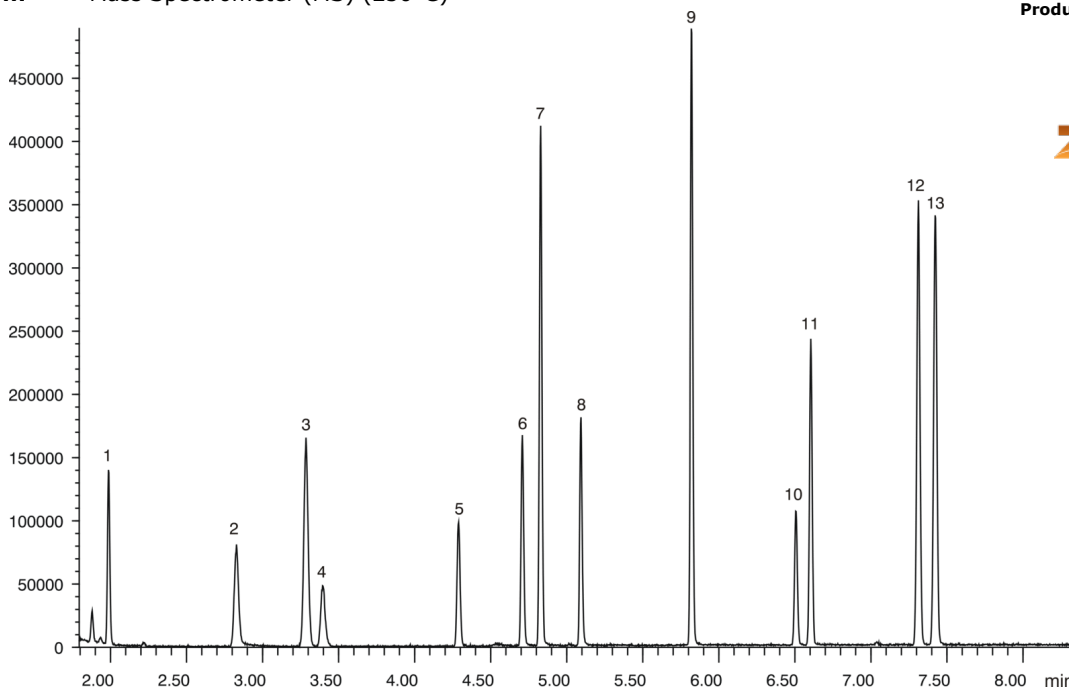