

Chlorthalidone

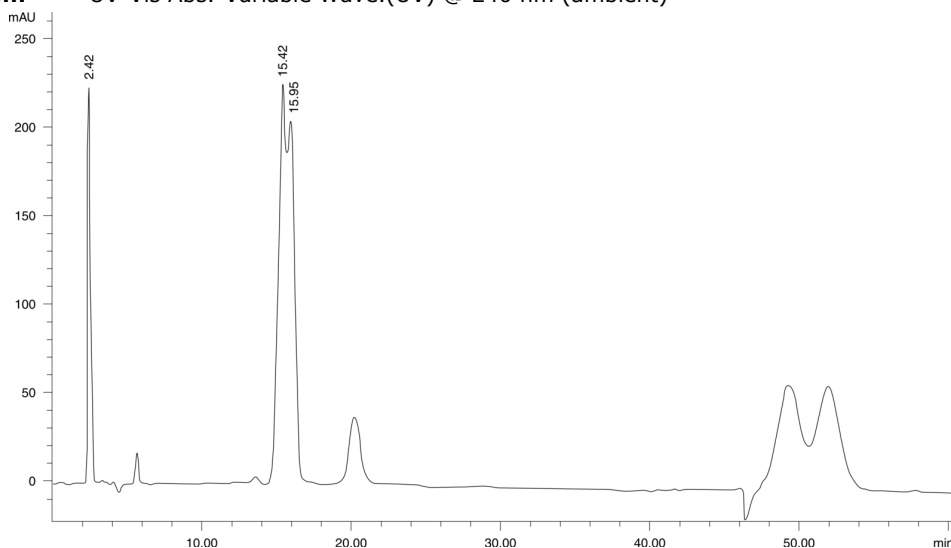
Column: Chirex® (S)-VAL and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3014-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (75:20:5) ethanol-TFA 20:1
Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 7 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 240 nm (ambient)



Products used in this application:

CHIREX®

ANALYTES:

- 1 Chlorthalidone
Retention Time: 15.42 min
- 2 Chlorthalidone
Retention Time: 15.95 min

