


## Dichloroisoproterenol

**Column:** Chirex® (S)-VAL and (R)-NEA, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-3014-E0  
**Elution Type:** Isocratic  
**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA 77:20:3 ethanol-TFA 20:1  
**Gradient**  


| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

  
**Flow Rate:** 0.7 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 272 nm (ambient)

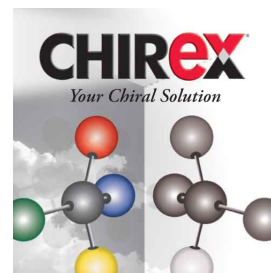

Products used in this application:

**CHIREX**®

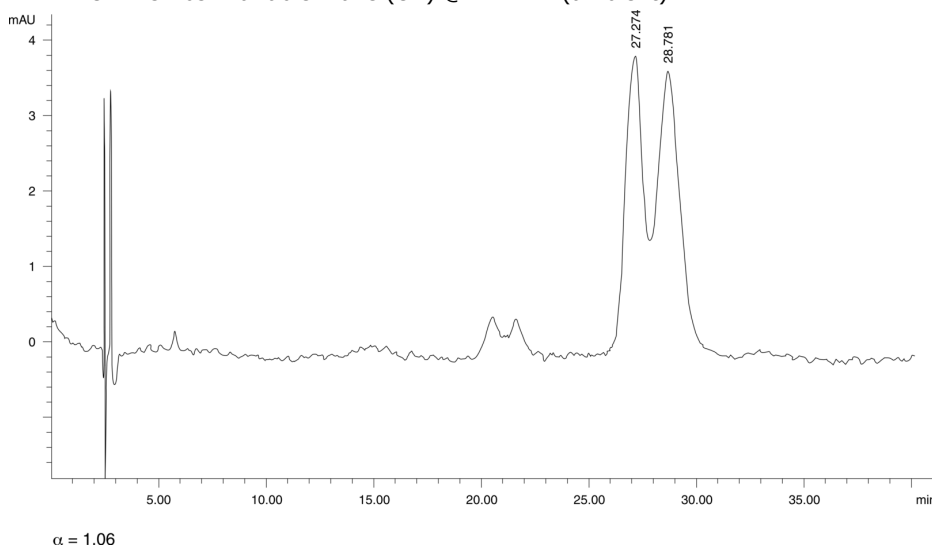

### ANALYTES:

- 1 Dichloroisoproterenol  
Retention Time: 27.274 min
- 2 Dichloroisoproterenol  
Retention Time: 28.781 min

