

## Isoetharine

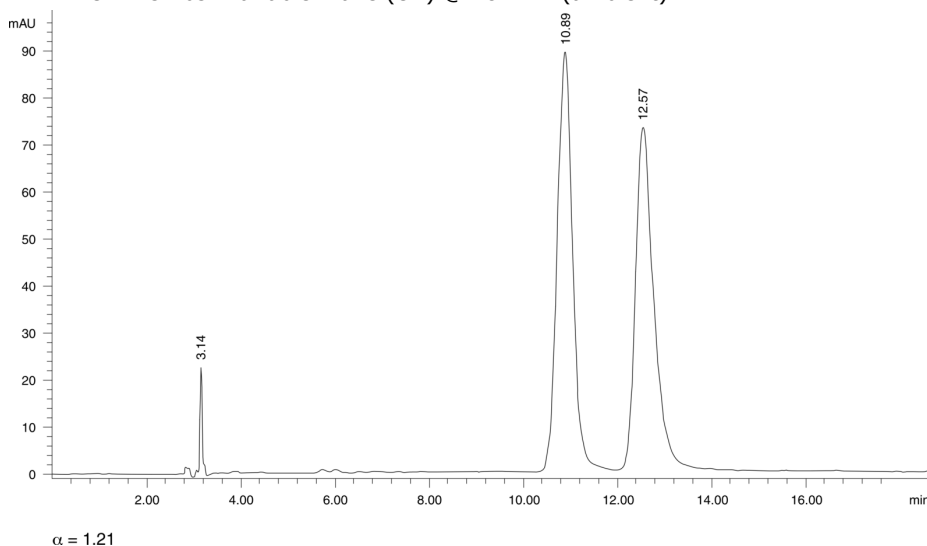
**Column:** Chirex® (S)-LEU and (R)-NEA, LC Column 250 x 4.6 mm, Ea  
**Dimensions:** 250 x 4.6 mm ID  
**Order No:** 00G-3020-E0  
**Elution Type:** Isocratic  
**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (55:35:10) ethanol-TFA 20:1  
**Gradient**  

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

  
**Flow Rate:** 1 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 282 nm (ambient)



Products used in this application:

**CHIREX®**

### ANALYTES:

- 1 Isoetharine  
Retention Time: 10.89 min
- 2 Isoetharine  
Retention Time: 12.57 min

