

Labetalol

Column: Chirex® (S)-ICA and (R)-NEA, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-3022-E0

Elution Type: Isocratic


Eluent A: Hexane/Dichloroethane/Ethanol-TFA (55:35:10)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 308 nm (ambient)

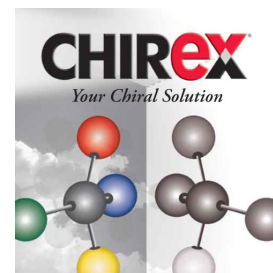

Products used in this application:

CHIREX®

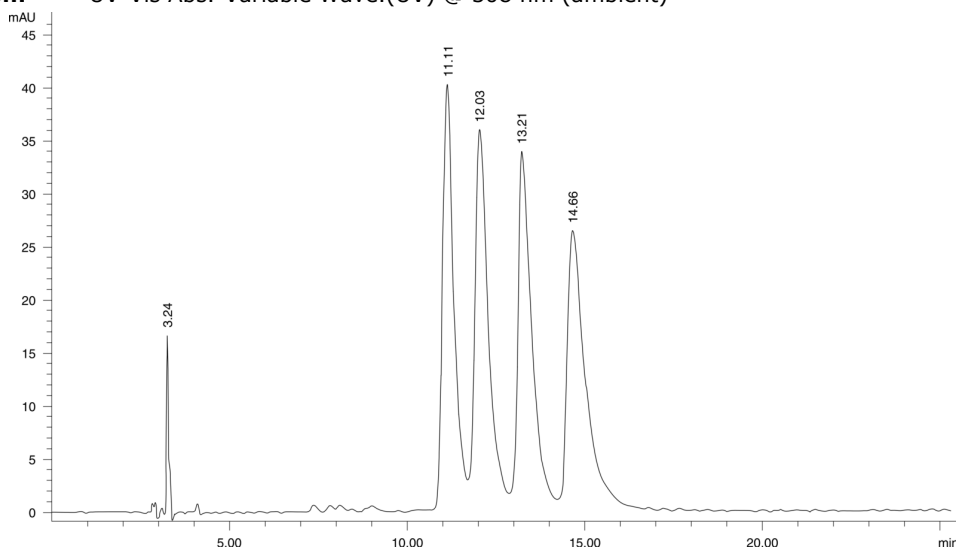

ANALYTES:

- 1 Labetalol
Retention Time: 11.11 min
- 2 Labetalol
Retention Time: 12.03 min
- 3 Labetalol
Retention Time: 13.21 min
- 4 Labetalol
Retention Time: 14.66 min

