

HPLC Application

ID No.: 13682

Primaquine

Column: Chirex[®] (S)-VAL and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3014-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (55:35:10) ethanol/TFA (20:1)
Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Refractive Index @ 0.0000000000 (ambient)



Products used in this application:

CHIREX[®]

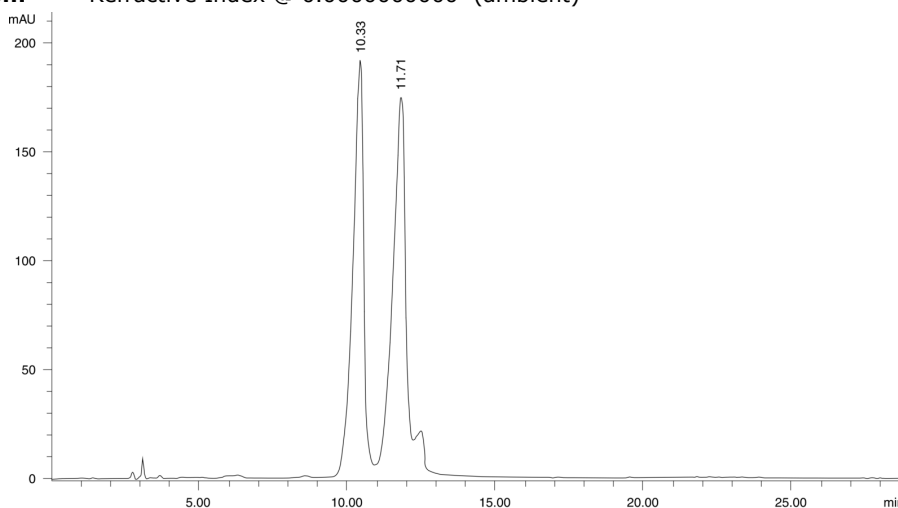

$\alpha = 1.18$

ANALYTES:

- 1 Primaquine
Retention Time: 10.33 min
- 2 Primaquine
Retention Time: 11.71 min

