

HPLC Application

ID No.: 14034

2-Amino-n-butyrac acid

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic


Eluent A: 1 mM Copper (II) sulfate in water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

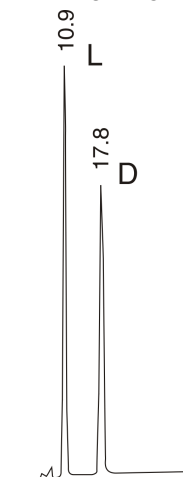

$$\alpha = 1.80$$

ANALYTES:

- 1** L-2-Amino-n-butyrac acid
Retention Time: 10.85 min
- 2** 2-Amino-n-butyrac acid
Retention Time: 17.703 min

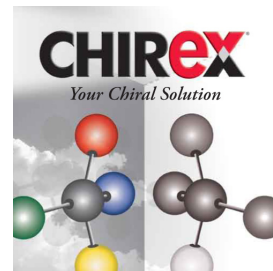

Products used in this application:

CHIREX®

