

EPA Method 554, DNPH derivatized Carbonyl Compounds

Column: Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4252-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	30	70
	2	20	30	70
	3	35	0	100
	4	45	0	100

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: Refractive Index @ 0.0000000000 (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

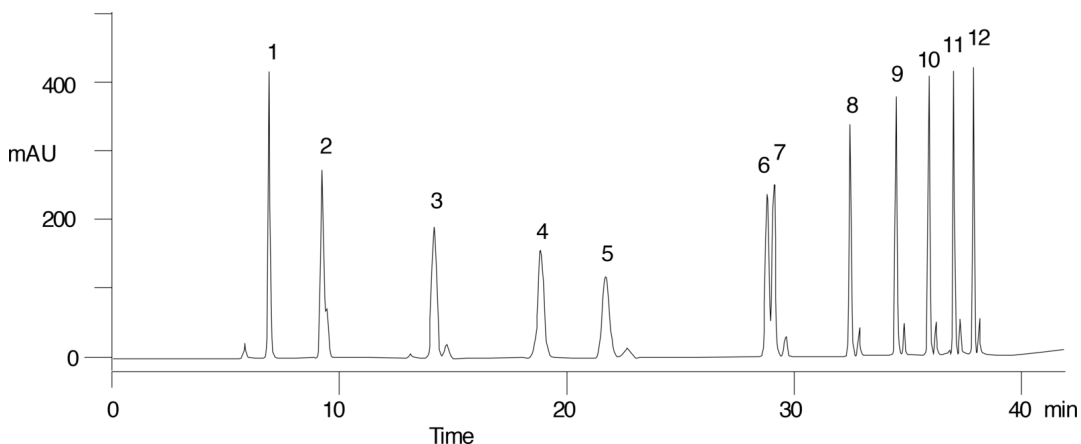

EPA Method 554, DNPH derivatized Carbonyl Compounds

ANALYTES:

- 1** Formaldehyde
- 2** Acetaldehyde
- 3** Propanal
- 4** Butanal
- 5** Pentanal
- 6** Hexanal
- 7** Heptanal
- 8** Octanal
- 9** Nonanal
- 10** Decanal
- 11** Cyclohexanone
- 12** Crotonaldehyde

