

## siRNA on Jupiter Proteo

**Column:** Jupiter® 4 µm Proteo 90 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4396-E0

**Elution Type:** Gradient

**Eluent A:** 100mM TEAA pH 7

**Eluent B:** ACN

| Gradient | Step No. | Time (min) | Pct A | Pct B |
|----------|----------|------------|-------|-------|
| Profile: | 1        | 0.5        | 96    | 4     |
|          | 2        | 25         | 88    | 12    |

**Flow Rate:** 1 mL/min


**Col. Temp.:** 30 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

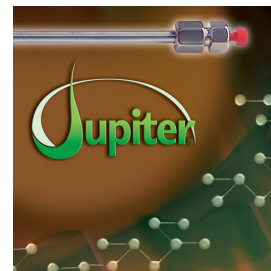

Products used in this application:

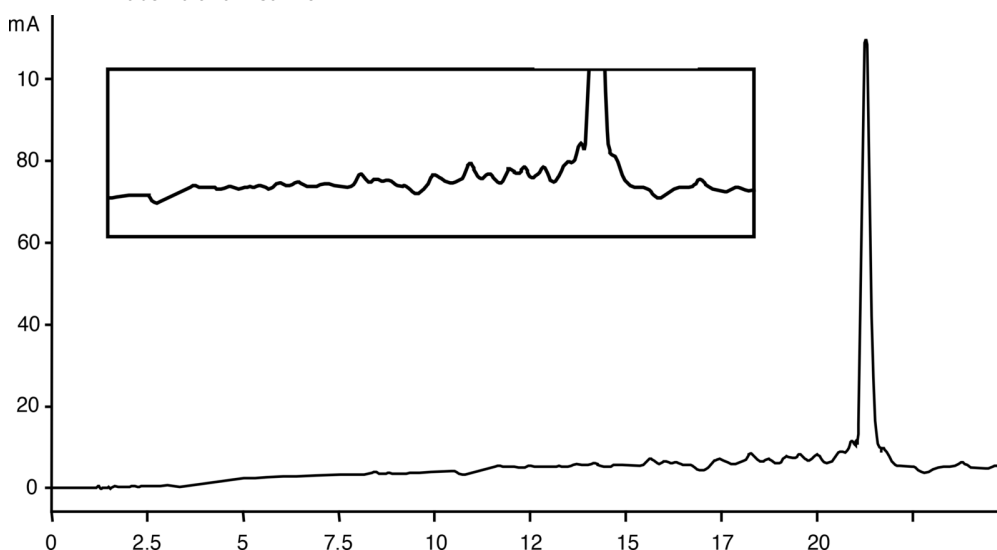

### ANALYTES:

1 siRNA 19mer

