


EPA Method 8100: PAH's on ZB-5 30x0.25x0.5

Column: Zebtron™ ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.50 µm, Ea
Phase: 5% Phenyl 95% Dimethylpolysiloxane
Dimensions: 30 meters x 0.25 mm x 0.5 µm
Order No: 7HG-G002-17
Oven Profile: 170°C to 240°C at 15°C/min to 310°C at 8°C/min to 330°C at 20°C/min (hold 5 min)
Carrier Gas: Constant Flow Helium, 1.6 mL/min
Injection: split 20:1 1 µL @ 260°C
Detection: Flame Ionization (FID) (325°C)

14778



Products used in this application:

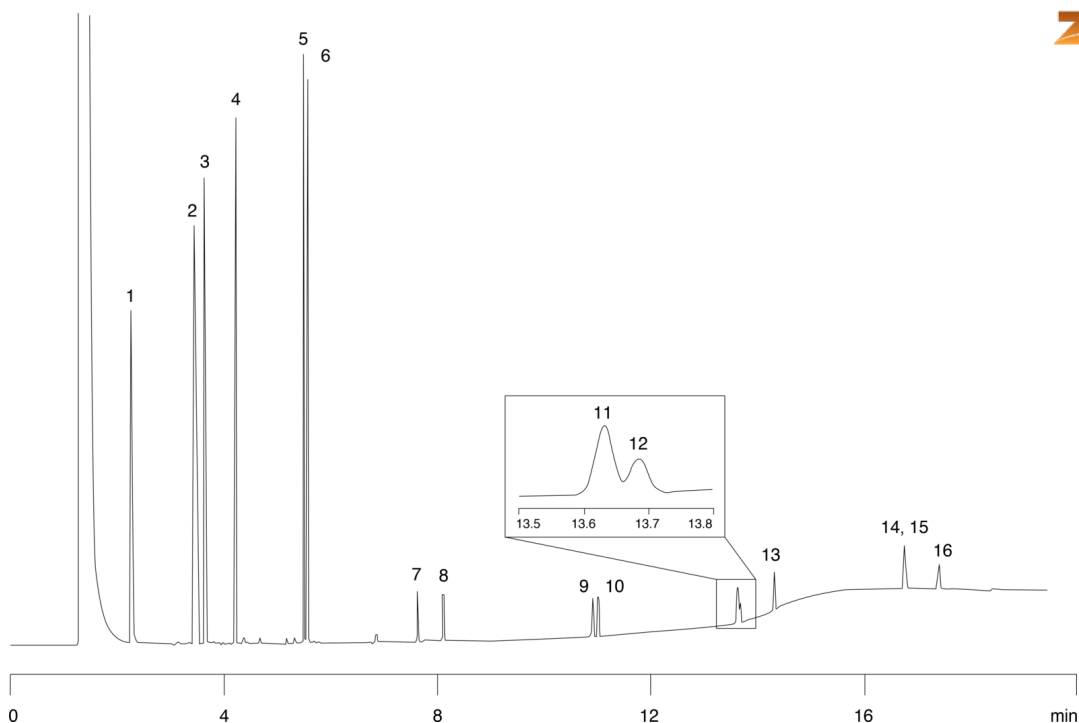

EPA Method 8100: PAH's on ZB-5 30x0.25x0.5

ANALYTES:

- 1 Naphthalene
- 2 Acenaphthylene
- 3 Acenaphthene
- 4 Fluorene
- 5 Phenanthrene
- 6 Anthracene
- 7 Fluoranthene
- 8 Pyrene
- 9 Benz[a]anthracene
- 10 Chrysene
- 11 Benzo[b]fluoranthene
- 12 Benzo[k]fluoranthene
- 13 Benzo[a]pyrene
- 14 Benzo[g,h,i]perylene
- 15 Dibenzo[a,h]anthracene
- 16 Indeno[1,2,3-cd]pyrene

