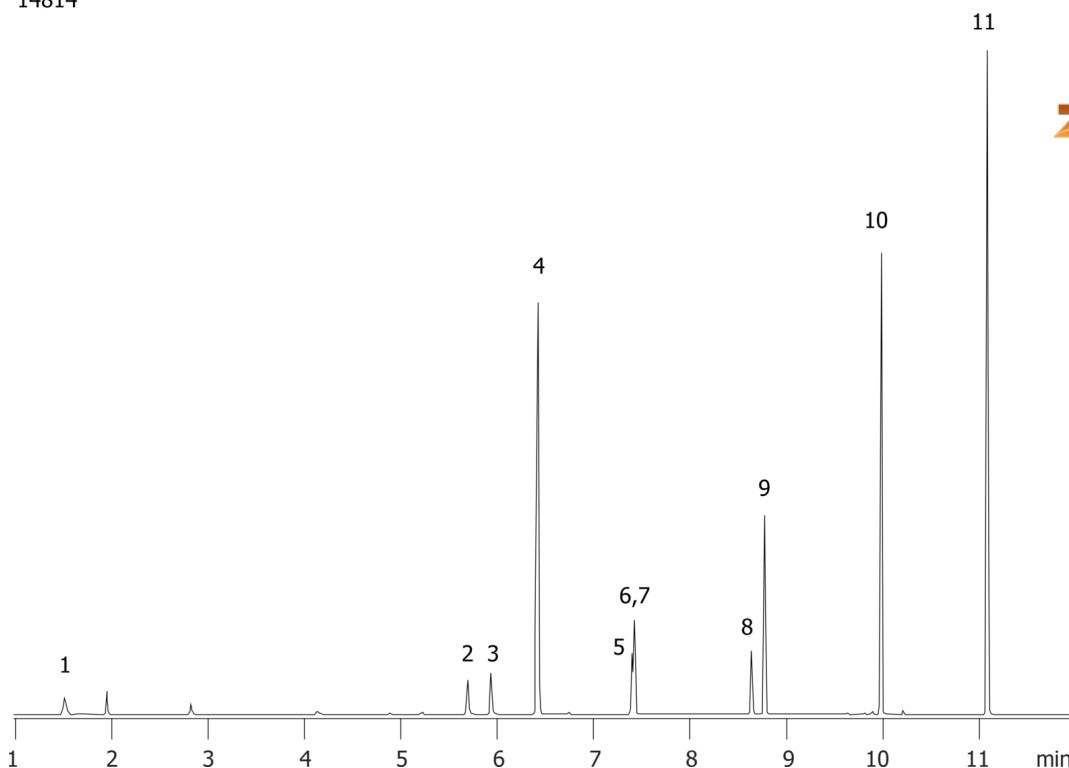


EPA Method 552.2: Haloacetic Acids on ZB-5 30x0.25x0.25 (fast)**Column:** Zebtron[™] ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G002-11**Oven Profile:** 40°C for 4 to 200°C at 14°C/min.**Carrier Gas:** Constant Flow Helium, 2.1 mL/min**Injection:** split 20:1 1 μ L @ 200°C**Detection:** Electron Capture (ECD) (260°C)

14814

**Products used in this application:**

EPA Method 552.2: Haloacetic Acids on ZB-5 30x0.25x0.25 (fast)

ANALYTES:

- 1 Chloroacetic acid
- 2 Bromoacetic acid
- 3 Dichloroacetic acid
- 4 Dalapon
- 5 1,2,3-Trichloropropane
- 6 Bromochloroacetic acid
- 7 Trichloroacetic acid
- 8 Dibromoacetic acid
- 9 Dibromodichloroacetic acid
- 10 Chlorodibromoacetic acid
- 11 Tribromoacetic acid

