


Nitrosamines ZB-5 30x0.25x0.5**Column:** Zebtron[™] ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.50 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.5 µm**Order No:** 7HG-G002-17**Oven Profile:** 40°C to 140°C at 6°C/min to 300°C at 25°C/min.**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** Split 50:1 1 µL @ 200°C**Detection:** Mass Selective (MSD) (240°C)**Analyst Note:** All analytes are 200 µg/mL**Products used in this application:****ANALYTES:**

- 1 N-Nitrosodimethylamine
- 2 N-Nitrosomethylethylamine
- 3 N-Nitrosodiethylamine
- 4 N-Nitrosopyrrolidine
- 5 N-Nitrosomorpholine
- 6 N-Nitrosodi-n-propylamine
- 7 N-Nitrosopiperidine
- 8 N-Nitrosodi-n-butylamine
- 9 N-Nitrosodiphenylamine

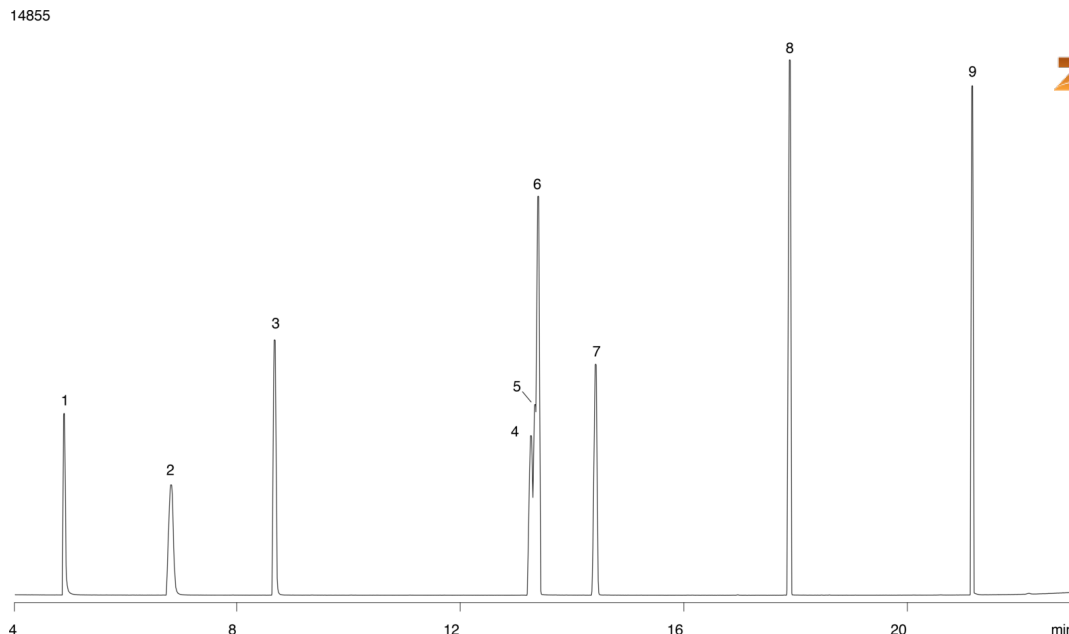