


**Unsaturated Fatty Acid Methyl Esters (FAME's) - ZB-FFAP 60x0.25x0.25****Column:** Zebron™ ZB-FFAP, GC Cap. Column 60 m x 0.25 mm x 0.25 µm, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 60 meters x 0.25 mm x 0.25 µm**Order No:** 7KG-G009-11**Oven Profile:** 200°C to 260°C at 2°C/min for 30 min.**Carrier Gas:** Constant Flow Helium, 2.4 mL/min**Injection:** Split 40:1 0.1 µL @ 220°C**Detection:** Refractive Index (250°C)**Analyst Note:** Recommended Accessories:  
Recommended Liner: Zebron PLUS Single Taper Z-Liner, 4 mm ID  
Liner Part No.: AG2-0A13-05 (for Agilent systems)**Products used in this application:**

14856

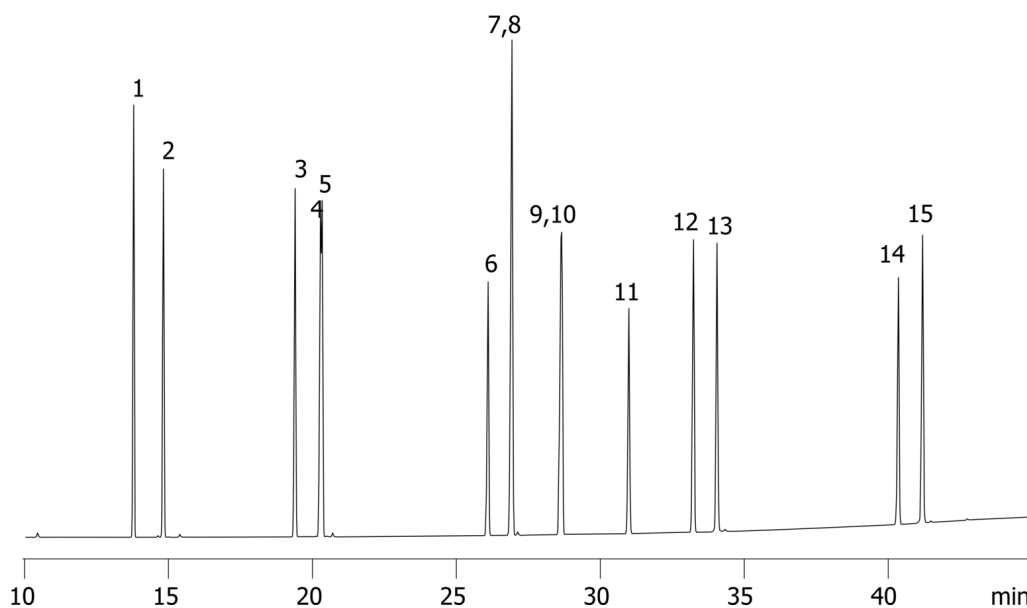

## Unsaturated Fatty Acid Methyl Esters (FAME's) - ZB-FFAP 60x0.25x0.25

### ANALYTES:

- 1 Methyl myristate
- 2 Methyl myristoleate
- 3 Methyl palmitate
- 4 Methyl palmitelaidate
- 5 Methyl palmitoleate
- 6 Methyl stearate
- 7 Methyl oleate
- 8 Methyl elaidate
- 9 Methyl linoleate
- 10 Methyl linelaidate
- 11 Methyl linolenate
- 12 Methyl arachidate
- 13 Methyl gondonate
- 14 Methyl behenate
- 15 Methyl erucate

