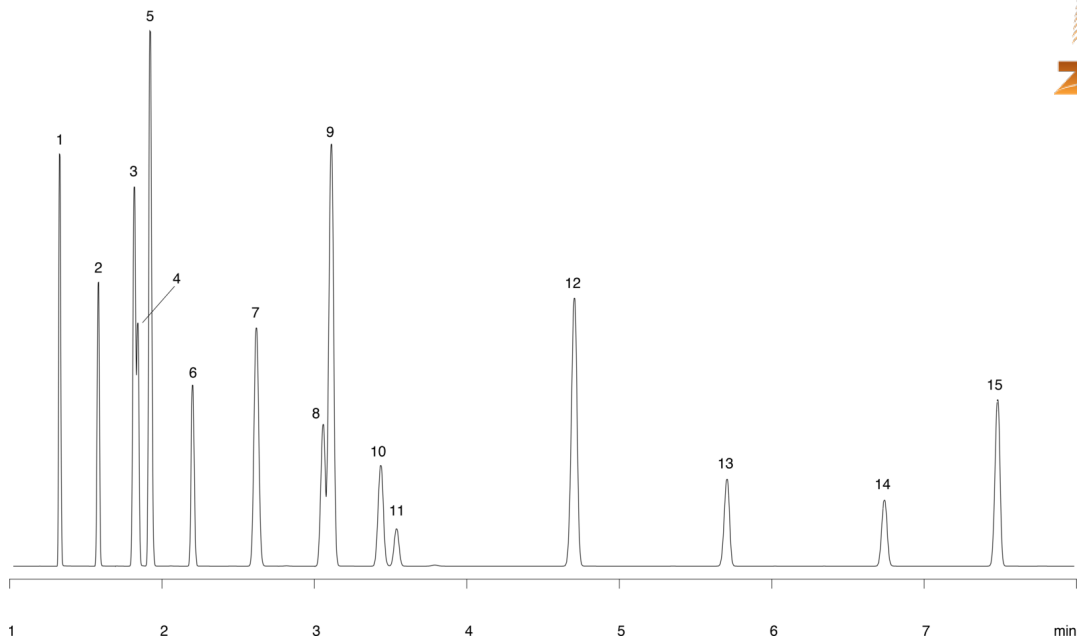


Abused Inhalants ZB-5 30x0.32x1.0**Column:** Zebron[™] ZB-5, GC Cap. Column 30 m x 0.32 mm x 1.00 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 1 µm**Order No:** 7HM-G002-22**Oven Profile:** 40°C for 4 min to 120°C at 10°C/min.**Carrier Gas:** Constant Flow Helium, 3.5 mL/min**Injection:** split 75:1 100 µL @ 200°C**Detection:** Flame Ionization (FID) (250°C)

14859

**Products used in this application:**

Abused Inhalants ZB-5 30x0.32x1.0

ANALYTES:

- 1** Methanol
- 2** Ethanol
- 3** Acetone
- 4** Isopropanol
- 5** Ether
- 6** Methylene chloride
- 7** Methyl tert-butyl ether
- 8** Methyl ethyl ketone
- 9** Hexane
- 10** Ethyl acetate
- 11** Chloroform
- 12** Benzene
- 13** Trichloroethylene
- 14** Methyl isobutyl ketone
- 15** Toluene

