


EPA Method 611 - Haloethers ZB-5 15x0.53x1.5**Column:** Zebron[™] ZB-5, GC Cap. Column 15 m x 0.53 mm x 1.50 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.53 mm x 1.5 μ m**Order No:** 7EK-G002-28**Oven Profile:** 100°C to 300°C at 15°C/min for 10 min.**Carrier Gas:** Constant Flow Helium, 3 mL/min**Injection:** split 20:1 0.2 μ L @ 200°C**Detection:** Flame Ionization (FID) (300°C)

14864

**Products used in this application:****ANALYTES:**

- 1 bis(2-Chloroethyl)ether
- 2 bis(2-Chloroisopropyl)ether
- 3 bis(2-Chloroethoxy)ether
- 4 4-Chlorophenylphenylether
- 5 4-Bromophenyl phenyl ether

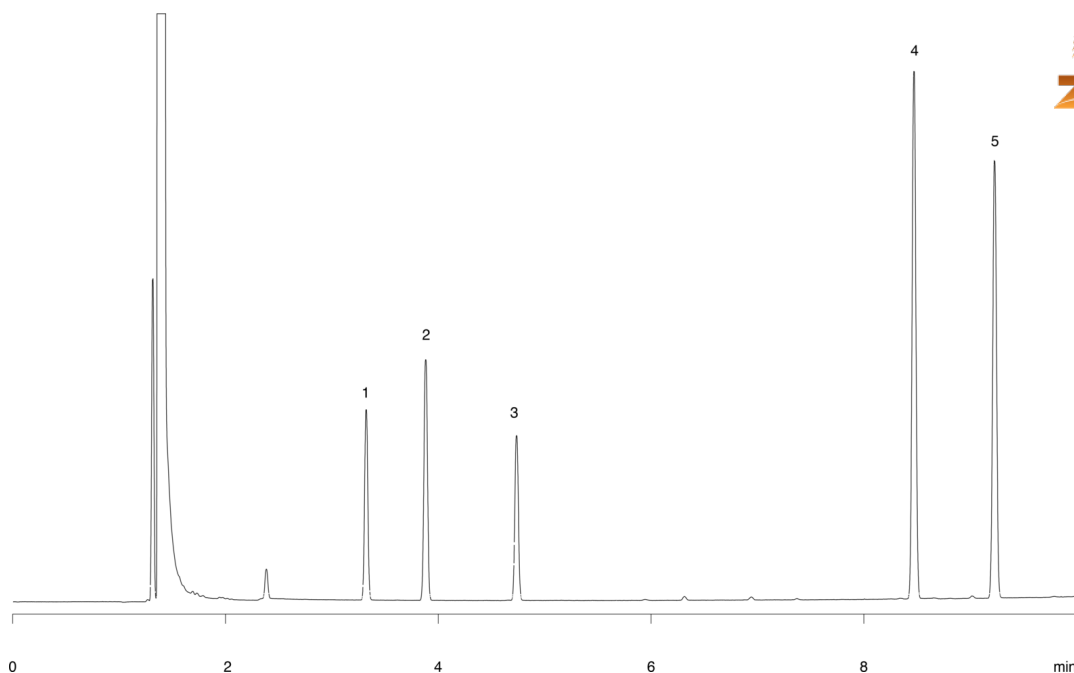