


**EPA Method 607 - Nitrosamines: ZB-Wax 30x0.32x0.5****Column:** Zebtron<sup>TM</sup> ZB-WAX, GC Cap. Column 30 m x 0.32 mm x 0.50 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.32 mm x 0.5 µm**Order No:** 7HM-G007-17**Oven Profile:** 40°C for 4 min to 240°C at 20°C/min for 10 min.**Carrier Gas:** Constant Flow Helium, 7 mL/min**Injection:** split 15:1 1 µL @ 240°C**Detection:** Refractive Index (240°C)

14869

**Products used in this application:****ANALYTES:**

- 1 N-Nitrosodimethylamine
- 2 N-Nitrosodi-n-propylamine
- 3 N-Nitrosodiphenylamine

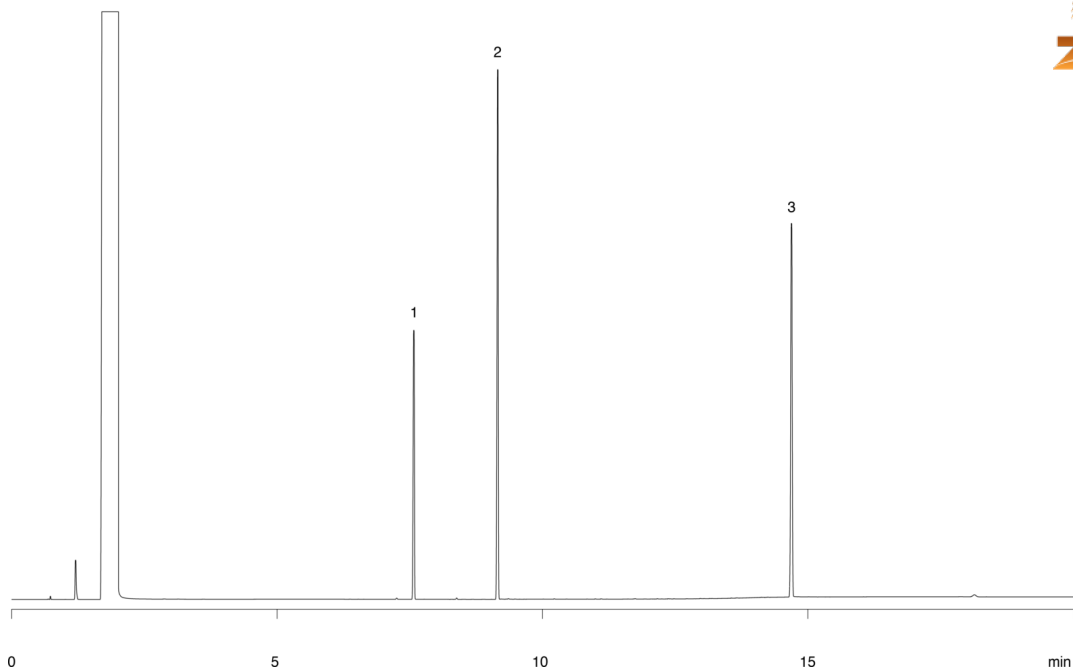