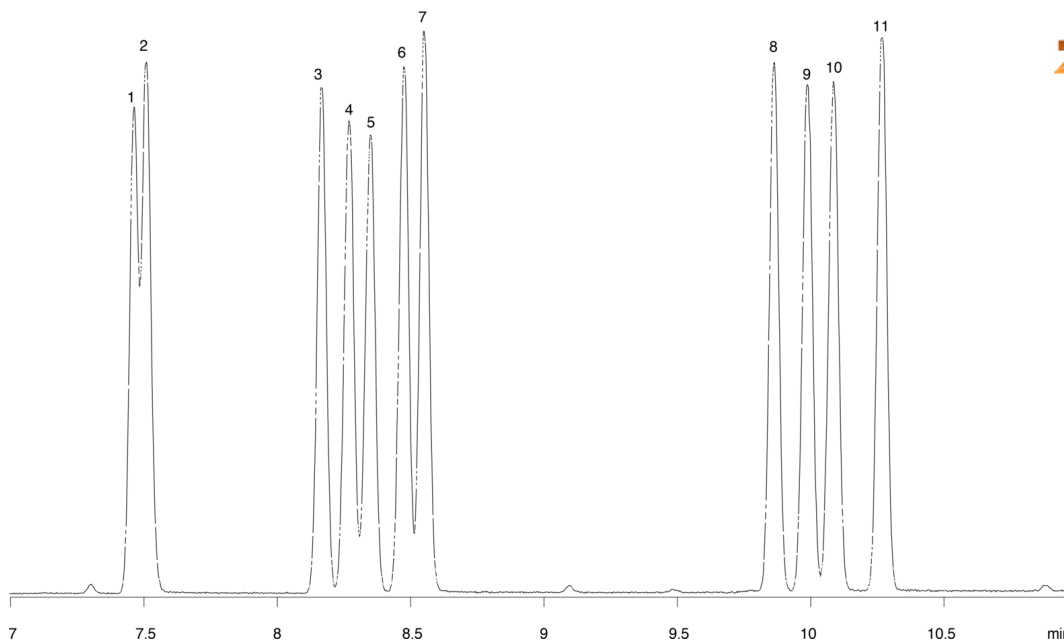


Triazine herbicides: EPA Method 619 ZB-1701 30x0.25x0.25**Column:** Zebtron™ ZB-1701, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 14% Cyanopropylphenyl 86% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G006-11**Oven Profile:** 170°C to 250°C at 6°C/min.**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** split 40:1 1 µL @ 250°C**Detection:** Mass Selective (MSD) (240°C)

14873

**Products used in this application:**

Triazine herbicides: EPA Method 619 ZB-1701 30x0.25x0.25

ANALYTES:

- 1** Prometon
- 2** Atraton
- 3** Propazine
- 4** Atrazine
- 5** Simazine
- 6** Terbutylazine
- 7** Secbumetone
- 8** Prometryn
- 9** Ametryn
- 10** Simetryn
- 11** Tertbutryn

