


Primary Amines - ZB-Wax 30x0.53x1.0**Column:** Zebron[™] ZB-WAX, GC Cap. Column 30 m x 0.53 mm x 1.00 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.53 mm x 1 μ m**Order No:** 7HK-G007-22**Oven Profile:** 40°C for 1 min to 150°C at 30°C/min for 2 min.**Carrier Gas:** Constant Flow Helium, 8 mL/min**Injection:** split 80:1 0.2 μ L @ 250°C**Detection:** Flame Ionization (FID) (250°C)

14877

**Products used in this application:****ANALYTES:**

- 1 Trimethylamine
- 2 Triethylamine
- 3 tert-Butyl amine
- 4 Diethylamine
- 5 Isobutylamine

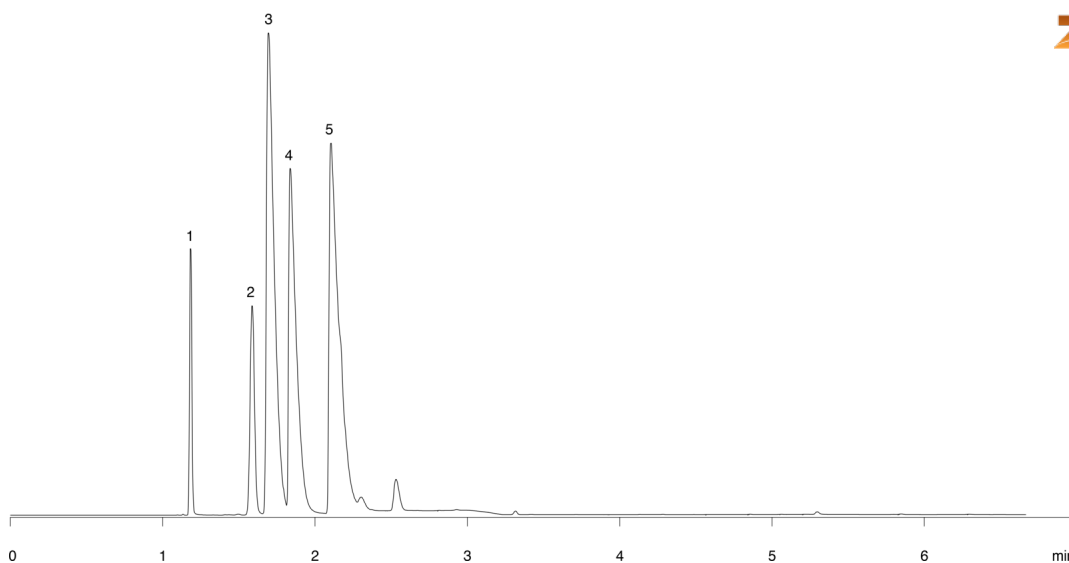