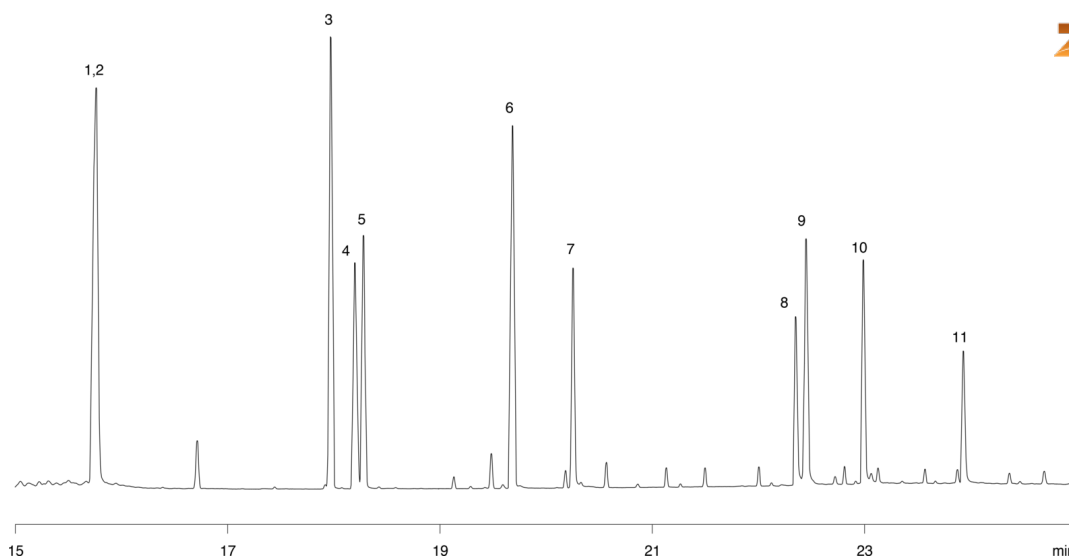


Phenols: EPA Method 604 - ZB-50 30x0.25x0.25**Column:** Zebtron[™] ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 50% Phenyl 50% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G004-11**Oven Profile:** 50°C for 10 min to 250°C at 15°C/min for 10 min.**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** split 20:1 1 μ L @ 280°C**Detection:** Flame Ionization (FID) (280°C)

14884

**Products used in this application:**

Phenols: EPA Method 604 - ZB-50 30x0.25x0.25

ANALYTES:

- 1** Phenol
- 2** 2-Chlorophenol
- 3** 2,4-Dimethylphenol
- 4** 2-Nitrophenol
- 5** 2,4-Dichlorophenol
- 6** 4-Chloro-3-methylphenol
- 7** 2,4,6-Trichlorophenol
- 8** 2,4-Dinitrophenol
- 9** 4-Nitrophenol
- 10** 2-Methyl-4,6-dinitrophenol
- 11** Pentachlorophenol

