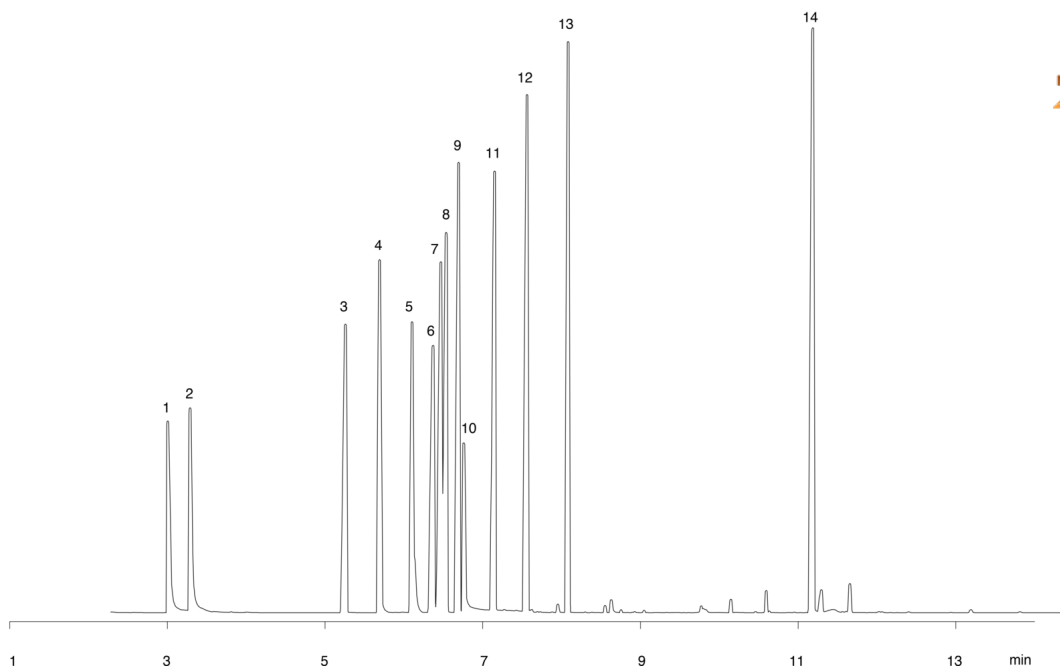


Aromatic Amines using 5 meter Guardian ZB-5 30x0.25x0.25**Column:** Zebron[™] ZB-5 w/GUARDIAN[™], 5 m, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G002-11-GGA**Oven Profile:** 60°C for 1min to 130°C at 12°C/min to 220°C at 30°C/min for 4 min.**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 15:1 1 μ L @ 220°C**Detection:** Refractive Index (240°C)

14944

**Products used in this application:**

Aromatic Amines using 5 meter Guardian ZB-5 30x0.25x0.25

ANALYTES:

- 1 Piperidine
- 2 2-Methylpiperidine
- 3 Aniline
- 4 Benzylamine
- 5 alpha-Phenylethylamine
- 6 N-Methylaniline
- 7 m-Toluidine
- 8 o-Toluidine
- 9 N,N-Dimethylaniline
- 10 beta-Phenylethylamine
- 11 N-Ethylaniline
- 12 2,4-Dimethylaniline
- 13 N,N-Diethylaniline
- 14 Dibenzylamine

