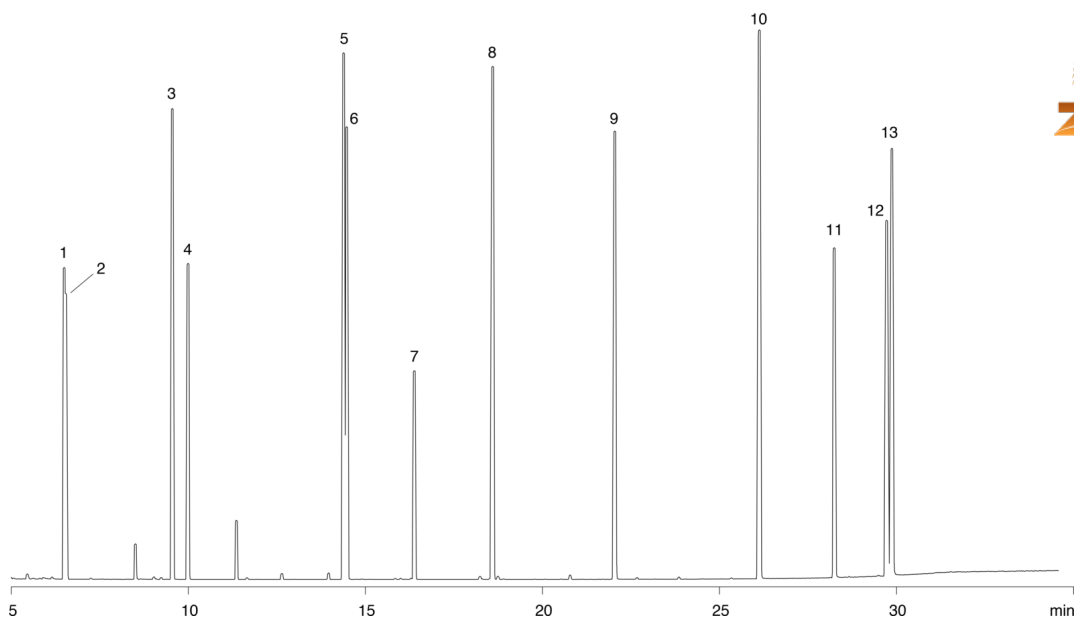


EPA Method 8270 Base Mix on ZB-5MS**Column:** Zebtron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:** 45°C to 320°C at 9°C/min for 4 min**Carrier Gas:** Constant Flow Helium, 1.24 mL/min**Injection:** Splitless :1 1 µL @ 250°C**Detection:** Mass Selective (MSD) (240°C)

14949

**Products used in this application:**

EPA Method 8270 Base Mix on ZB-5MS

ANALYTES:

- 1** 1,4-Dichlorobenzene-d4
- 2** 1,4-Dichlorobenzene
- 3** Naphthalene-d8
- 4** Hexachloro-1,3-butadiene
- 5** Acenaphthene-d10
- 6** Acenaphthene
- 7** N-Nitrosodiphenylamine
- 8** Phenanthrene-d10
- 9** Fluoranthene
- 10** Chrysene-d12
- 11** Di-n-octyl phthalate
- 12** Benzo[a]pyrene
- 13** Perylene-d12

