

Hydrocarbon Markers on ZB-5MS

Column: Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea

Phase: 5% Phenyl-Arylene 95% Dimethylpolysiloxane

Dimensions: 30 meters x 0.32 mm x 0.25 µm

Order No: 7HM-G010-11

Oven Profile: 40°C to 340°C at 5°C/min for 10 min

Carrier Gas: Constant Flow Helium, 2.2 mL/min


Injection: split 20:1 1 µL @ 360°C

Detection: Flame Ionization (FID) (360°C)

15424



Products used in this application:

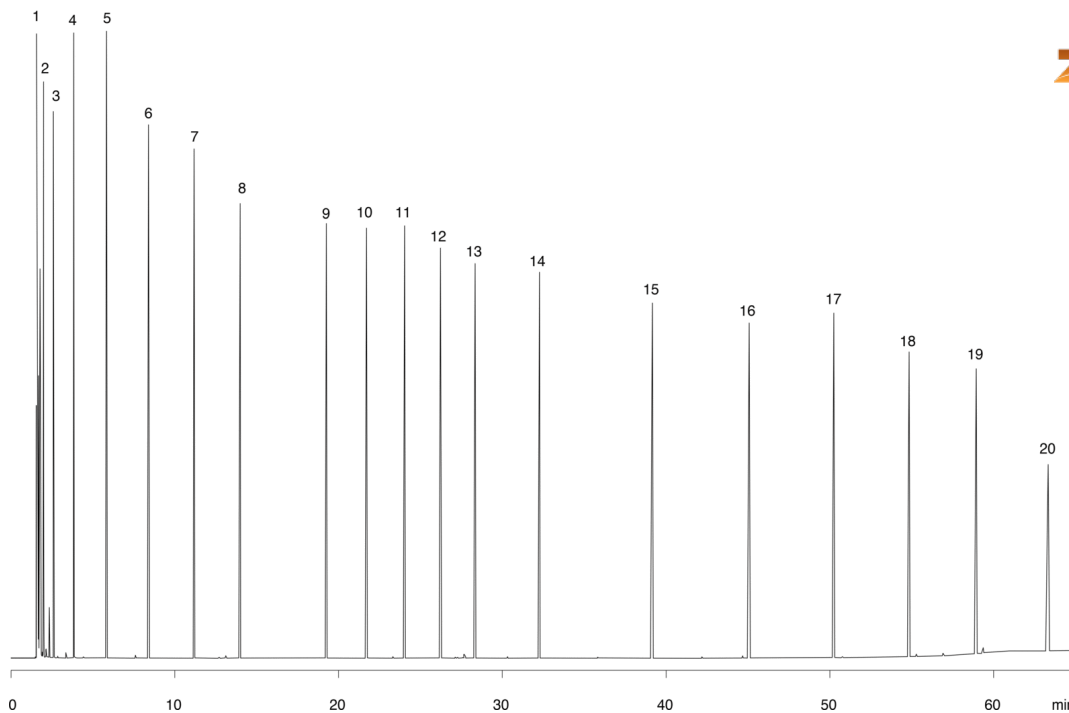

Hydrocarbon Markers on ZB-5MS

ANALYTES:

- 1 Pentane (C5)
- 2 Hexane (C6)
- 3 Heptane (C7)
- 4 Octane (C8)
- 5 Nonane (C9)
- 6 Decane (C10)
- 7 Undecane (C11)
- 8 Dodecane (C12)
- 9 Tetradecane (C14)
- 10 Pentadecane (C15)
- 11 Hexadecane (C16)
- 12 Heptadecane (C17)
- 13 Octadecane (C18)
- 14 Eicosane (C20)
- 15 Tetracosane (C24)
- 16 Octacosane (C28)
- 17 Dotriacontane (C32)
- 18 Hexatriacontane (C36)
- 19 Tetracontane (C40)
- 20 Tetratetracontane (C44)

