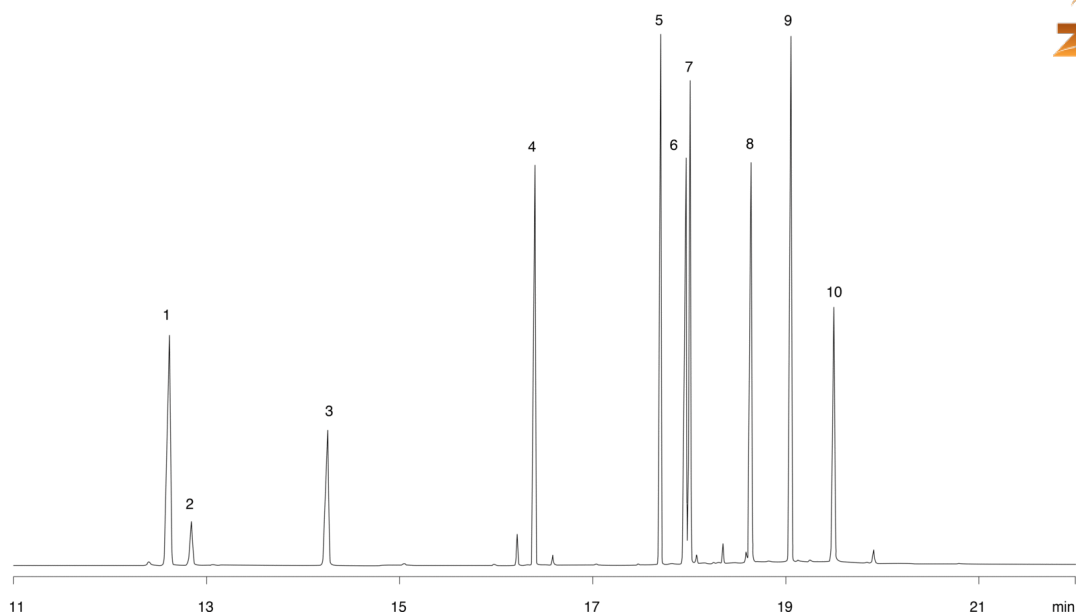


**Organo Phosphorous Pesticides Mix 1 on ZB-5MS**

**Column:** Zebron<sup>™</sup> ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea  
**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane  
**Dimensions:** 30 meters x 0.25 mm x 0.25 µm  
**Order No:** 7HG-G010-11  
**Oven Profile:** 100°C to 140°C at 20°C/min to 200°C at 5°C/min to 300°C at 25°C/min for 5 min  
**Carrier Gas:** Constant Flow Helium, 1 mL/min  
**Injection:** split 20:1 1 µL @ 250°C  
**Detection:** Mass Selective (MSD) (240°C)

15435

**Products used in this application:**

## Organo Phosphorous Pesticides Mix 1 on ZB-5MS

### ANALYTES:

- 1 Terbufos
- 2 Phosphamidon
- 3 Phosphamidon
- 4 Chlorfenvinphos
- 5 Ethion
- 6 Famphur
- 7 Carbophenothion
- 8 Phosmet
- 9 Leptophos
- 10 Azinphos-ethyl

