

Benzodiazepines (SPE-Plasma Blank) on Gemini C18**Column:** Gemini® 5 µm C18 110 Å, LC Column 150 x 3 mm, Ea**Dimensions:** 150 x 3 mm ID**Order No:** 00F-4435-Y0**Elution Type:** Gradient**Eluent A:** 0.1% Formic acid in Water**Eluent B:** 0.1% Formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	8	20	80
	3	9	20	80

Flow Rate: 0.8 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 2.0mm, 10/Pk Part No.: AJ0-7596

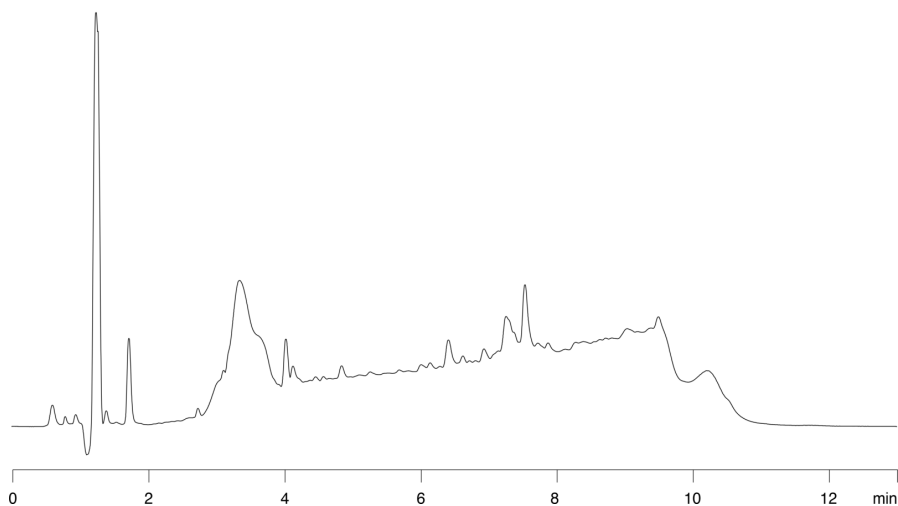
- Holder Part No.: KJ0-4282



Products used in this application:



15455

**ANALYTES:****1** blank