

Zearalone analysis using Synergi Hydro-RP

Column: Synergi™ 4 µm Hydro-RP 80 Å, LC Column 150 x 3 mm, Ea**Dimensions:** 150 x 3 mm ID**Order No:** 00F-4375-Y0**Elution Type:** Isocratic**Eluent A:** 0.1% Formic acid in Water**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50


Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, AQ C18 4 x 2.0mm, 10/Pk Part No.: AJ0-4286

- Holder Part No.: KJ0-4282



15531

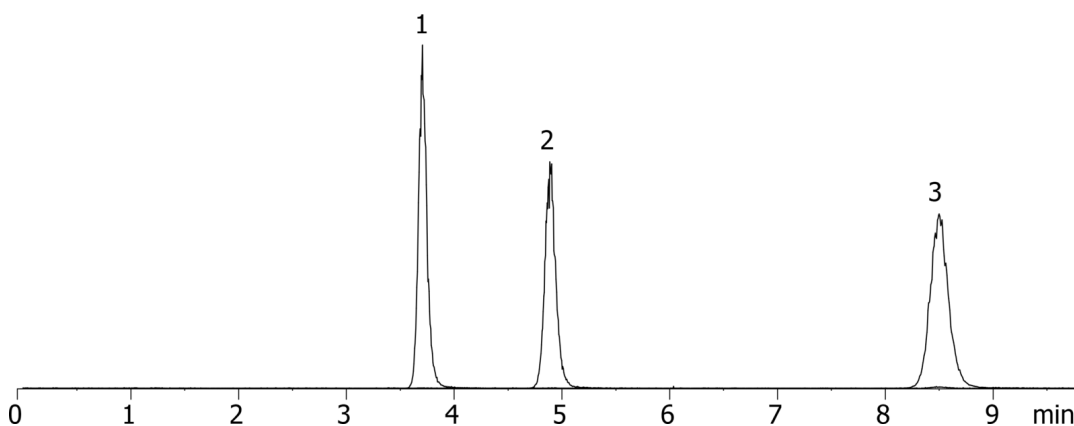

ANALYTES:

- 1 beta-Zearalenol
- 2 alpha-Zearalenol
- 3 Zearalenone

