# **HPLC Application** ID No.: **15623**



## Acid, Base, and Neutral Compounds by Luna 3μ C18(2) MercuryMS Cartridge 20 x 4.0mm

Luna® 3 µm C18(2) 100 Å, LC Column 20 x 4 mm, Mercury, Ea

20 x 4 mm ID **Dimensions:** Order No: 00M-4251-D0-CE

**Elution Type:** Gradient

Eluent A:

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	95
	2	1	5
	3	1.2	5
	4	1.21	95
	5	1 5	95

Flow Rate: 2.5 mL/min Col. Temp.: ambient

**Detection:** Refractive Index @ 0.000000000 (25 °C) **Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287


- Holder Part No.: KJ0-4282





Products used in this application:



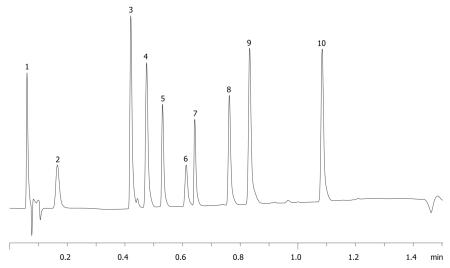

# **HPLC Application** ID No.: **15623**



## Acid, Base, and Neutral Compounds by Luna 3μ C18(2) MercuryMS Cartridge 20 x 4.0mm

### **ANALYTES:**

- Pyridine 1
- 2 Sulfanilamide
- 3 Quinidine
- Sulfathiazole 4
- 5 Triprolidine
- Phenol 6
- 7 Nortriptyline
- 3-Methyl-4-nitrobenzoic acid 8
- 9 5-Methyl-salicylaldehyde
- 10 Hexaphenone