


Drug Screen on Zebron ZB-1MS 20m x 0.18mm x 0.18um

Column: Zebron™ ZB-1MS, GC Cap. Column 20 m x 0.18 mm x 0.18 µm, Ea
Phase: 100% Dimethylpolysiloxane
Dimensions: 20 meters x 0.18 mm x 0.18 µm
Order No: 7FD-G011-08
Oven Profile: 50°C for 1 min to 125°C at 25°C/min to 200°C at 12°C/min to 325°C at 20°C/min for 0.5min
Carrier Gas: Constant Flow Helium, 1 mL/min
Injection: Split 20:1 1 µL @ 230°C
Detection: Flame Ionization (FID) (330°C)

15867



Products used in this application:

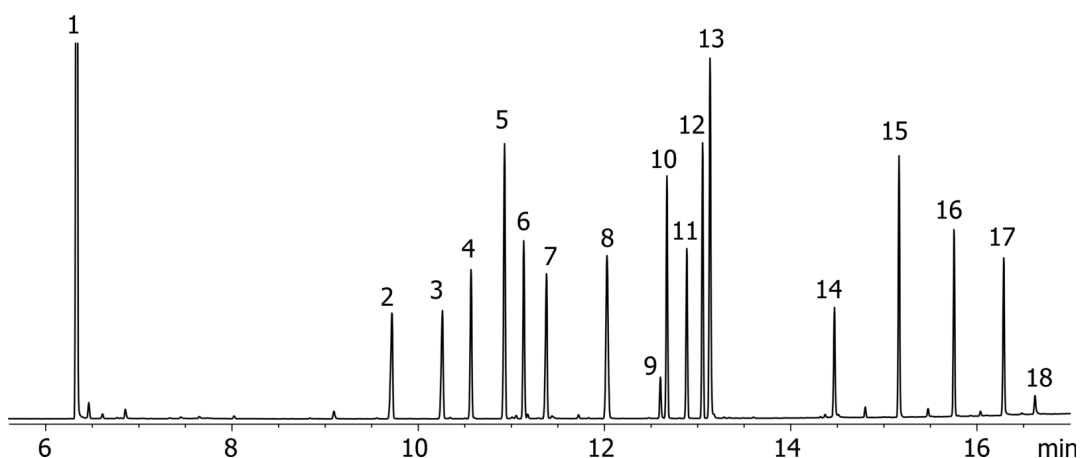

Drug Screen on Zebron ZB-1MS 20m x 0.18mm x 0.18um

ANALYTES:

- 1 Nicotine
- 2 Pentobarbital
- 3 Caffeine
- 4 Glutethimide
- 5 Lidocaine
- 6 Phencyclidine
- 7 Phenobarbital
- 8 EDDP
- 9 Methaqualone
- 10 Methadone
- 11 Cocaine
- 12 Trimipramine
- 13 Desipramine
- 14 Heroin
- 15 Fentanyl
- 16 Ibogaine
- 17 Triazolam
- 18 LSD

