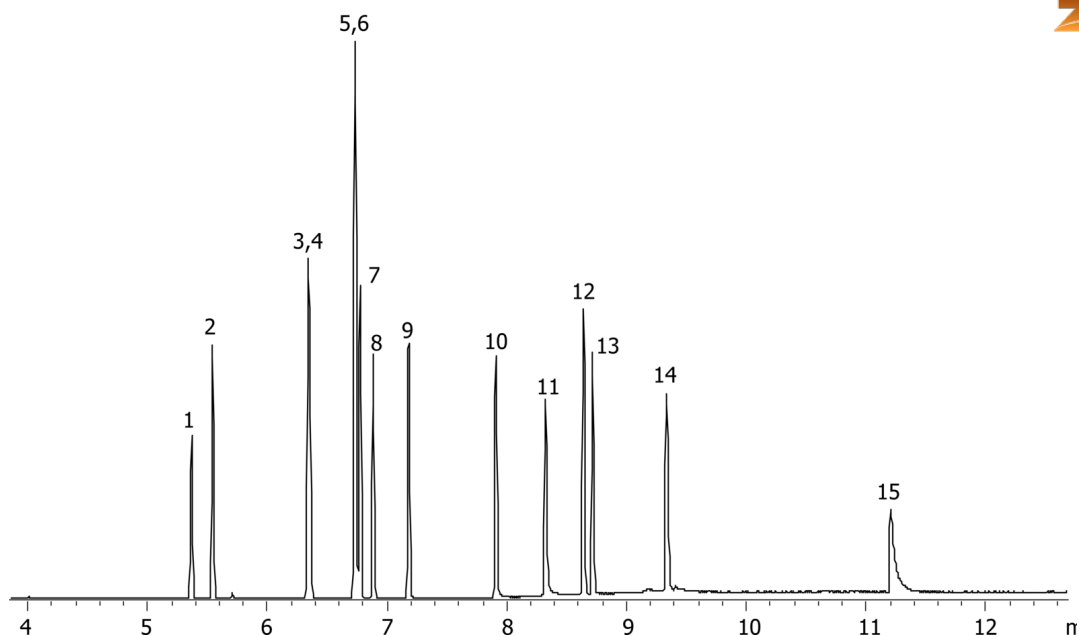


Phenols on Zebron ZB-WAX-Plus 10m x 0.10mm x 0.10µm**Column:** Zebron[™] ZB-WAXplus[™], GC Cap. Column 10 m x 0.10 mm x 0.10 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 10 meters x 0.1 mm x 0.1 µm**Order No:** 7CB-G013-02**Oven Profile:** 60°C to 240°C at 20°C/min for 6 min**Carrier Gas:** Constant Flow Helium, 0.23 mL/min**Injection:** Split 15:1 1 µL @ 220°C**Detection:** Mass Selective (MSD) (0°C)

15871

**Products used in this application:**

Phenols on Zebron ZB-WAX-Plus 10m x 0.10mm x 0.10um

ANALYTES:

- 1** 2-Nitrophenol
- 2** 2-Chlorophenol
- 3** Phenol
- 4** 2-Methylphenol
- 5** 4-Methylphenol
- 6** 2,4-Dimethylphenol
- 7** 3-Methylphenol
- 8** 2,6-Dichlorophenol
- 9** 2,4-Dichlorophenol
- 10** 2,4,5-Trichlorophenol
- 11** Benzoic acid
- 12** 4-Chloro-3-methylphenol
- 13** 2,4,6-Trichlorophenol
- 14** 2,3,4,6-Tetrachlorophenol
- 15** Pentachlorophenol

