

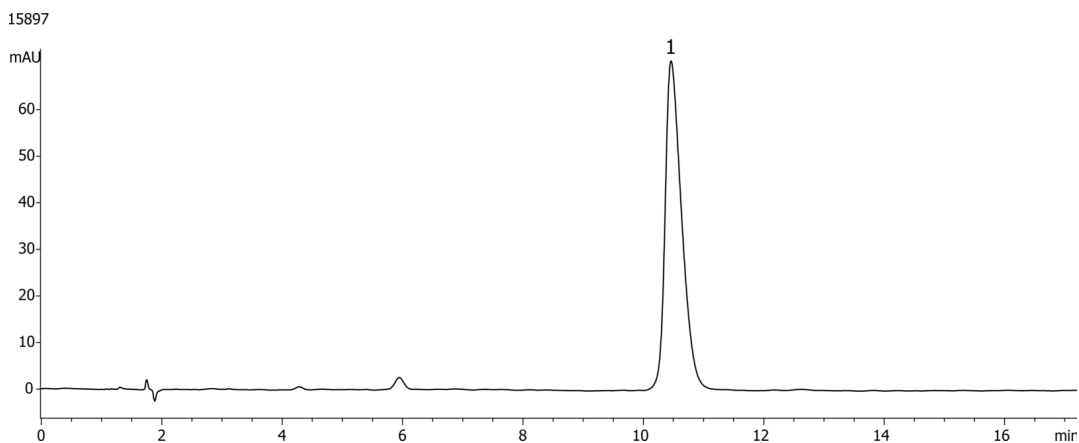
Oxacillin Alternate USP Method on Gemini 5µm C6 Phenyl**Column:** Gemini® 5 µm C6-Phenyl 110 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4444-E0**Elution Type:** Isocratic**Eluent A:** 0.1% Formic Acid in Water**Eluent B:** Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	35	65

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C6-Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-7915

- Holder Part No.: KJ0-4282

**Products used in this application:****ANALYTES:****1** Oxacillin