


Common Solvents on ZB-1MS**Column:** Zebtron[™] ZB-1MS, GC Cap. Column 30 m x 0.25 mm x 1.00 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 1 μ m**Order No:** 7HG-G011-22**Oven Profile:** 50°C for 4.5 min to 120°C at 20°C/min**Carrier Gas:** Constant Flow Helium, 3.2 mL/min**Injection:** Split 15:1 200 μ L @ 120°C**Detection:** Refractive Index (250°C)**Analyst Note:** analyte concentration was 0.05mg of each analyte per 20mL headspace vial**Products used in this application:**

Common Solvents on ZB-1MS

ANALYTES:

- 1 Methanol
- 2 Ethanol
- 3 Acetone
- 4 Isopropyl alcohol
- 5 n-Propanol
- 6 Methyl ethyl ketone
- 7 Ethyl acetate
- 8 Isopropyl acetate
- 9 Methoxypropanol (Propylene glycol monomethyl ether)
- 10 n-Propyl acetate
- 11 Methyl isobutyl ketone
- 12 3-Ethoxy propanol (Propylene glycol β -monoethyl ether)
- 13 Toluene
- 14 2,4-Pentanedione
- 15 n-Butyl acetate

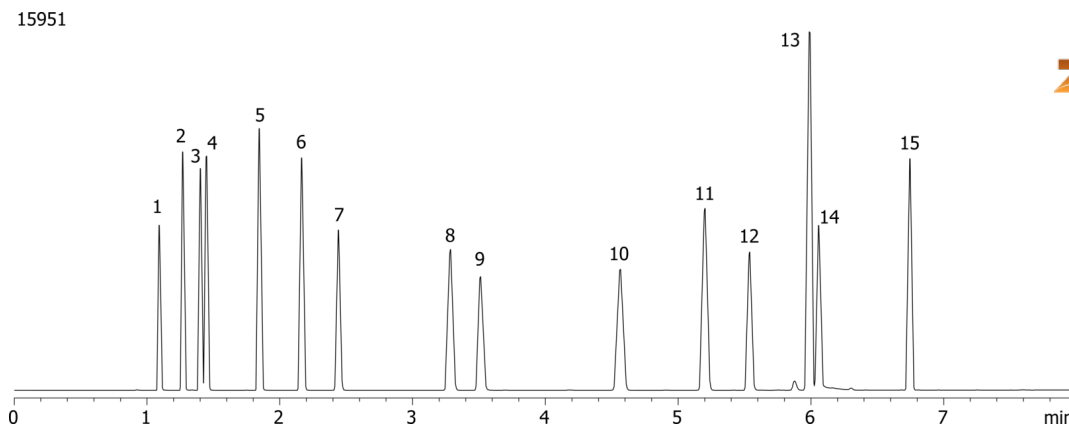