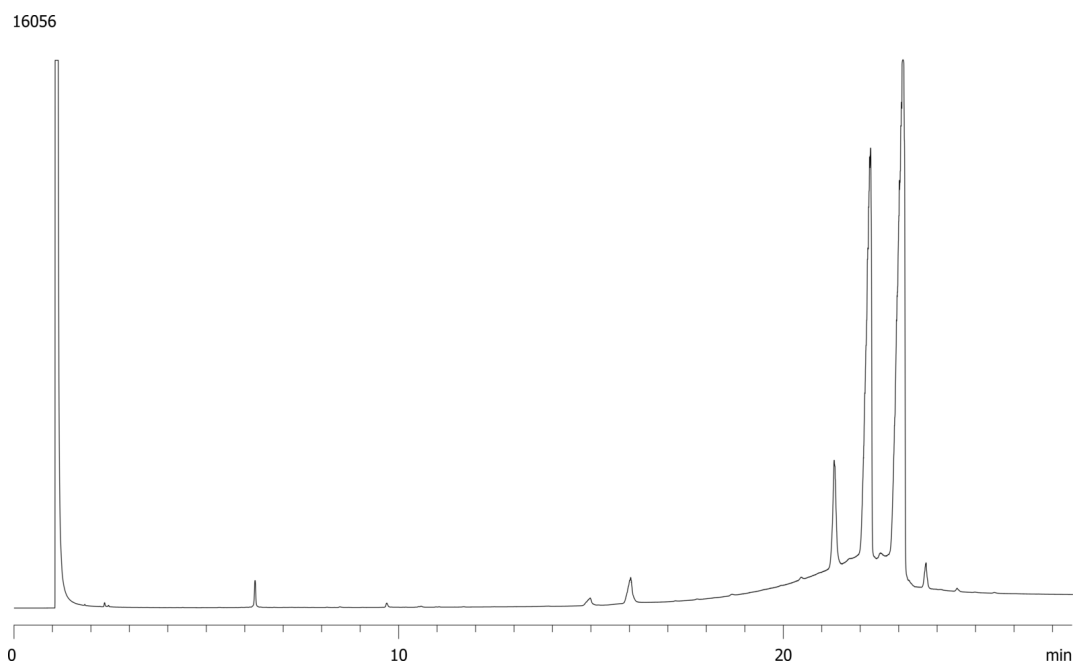


Olive Oil on ZB-5HT 30x0.25x0.10**Column:** Zebtron[™] ZB-5HT, GC Cap. Column 30 m x 0.25 mm x 0.10 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.1 μ m**Order No:** 7HG-G015-02**Oven Profile:** 220°C for 1 min to 400°C at 8°C/min for 4 min**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** On-Column :1 0.1 μ L @ 223°C**Detection:** Flame Ionization (FID) (415°C)**Analyst Note:** Sample was 1% in dichloromethane.
Recommended Accessories:
Recommended Liner: Zebtron PLUS Liner for Agilent & Thermo
4mm ID Direct Connect

Liner Part No.: AG2-0A50-05

**Products used in this application:****ANALYTES:****1** Olive oil