

Steroids on Gemini 5µm C18 250 x 3mm ID

Column: Gemini® 5 µm C18 110 Å, LC Column 250 x 3 mm, Ea

Dimensions: 250 x 3 mm ID

Order No: 00G-4435-Y0

Elution Type: Gradient

Eluent A: 0.1% Formic Acid in Water

Eluent B: 0.1% Formic Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40
	2	4	60	40
	3	15	25	75
	4	15.1	60	40
	5	20	60	40

Flow Rate: 0.6 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C18 4 x 2.0mm, 10/Pk Part No.: AJ0-7596

- Holder Part No.: KJ0-4282

16066

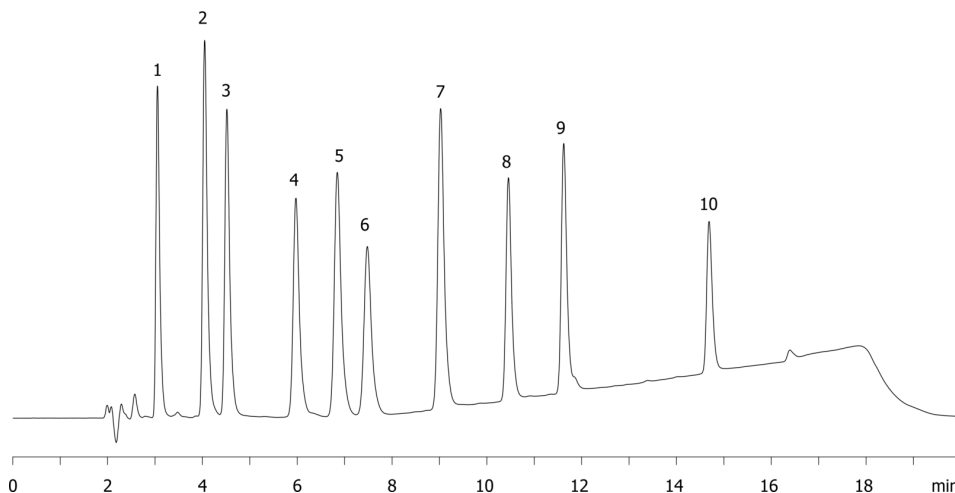

Products used in this application:



Steroids on Gemini 5µm C18 250 x 3mm ID

ANALYTES:

- 1** Triamcinolone
- 2** Prednisolone
- 3** Cortisone
- 4** Betamethasone
- 5** Corticosterone
- 6** Triamcinolone acetonide
- 7** 11-alpha-Hydroxyprogesterone
- 8** Cortisone acetate
- 9** 11-Ketoprogesterone
- 10** Betamethasone Valerate

