


Blood Alcohols on ZB-624**Column:** Zebron[™] ZB-624, GC Cap. Column 30 m x 0.32 mm x 1.80 µm, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 1.8 µm**Order No:** 7HM-G005-31**Oven Profile:** 40°C Isothermal**Carrier Gas:** Constant Flow Helium, 7 mL/min**Injection:** Splitless 13:1 200 µL @ 200°C**Detection:** Flame Ionization (FID) (220°C)**Analyst Note:** Static Headspace: 1mL sample at 50ppm in 20mL vial at 70°C for 15min.**Products used in this application:****ANALYTES:**

- 1 Acetaldehyde
- 2 Methanol
- 3 Ethanol
- 4 Acetone
- 5 Isopropanol
- 6 t-Butanol
- 7 n-Propanol
- 8 2-Butanol

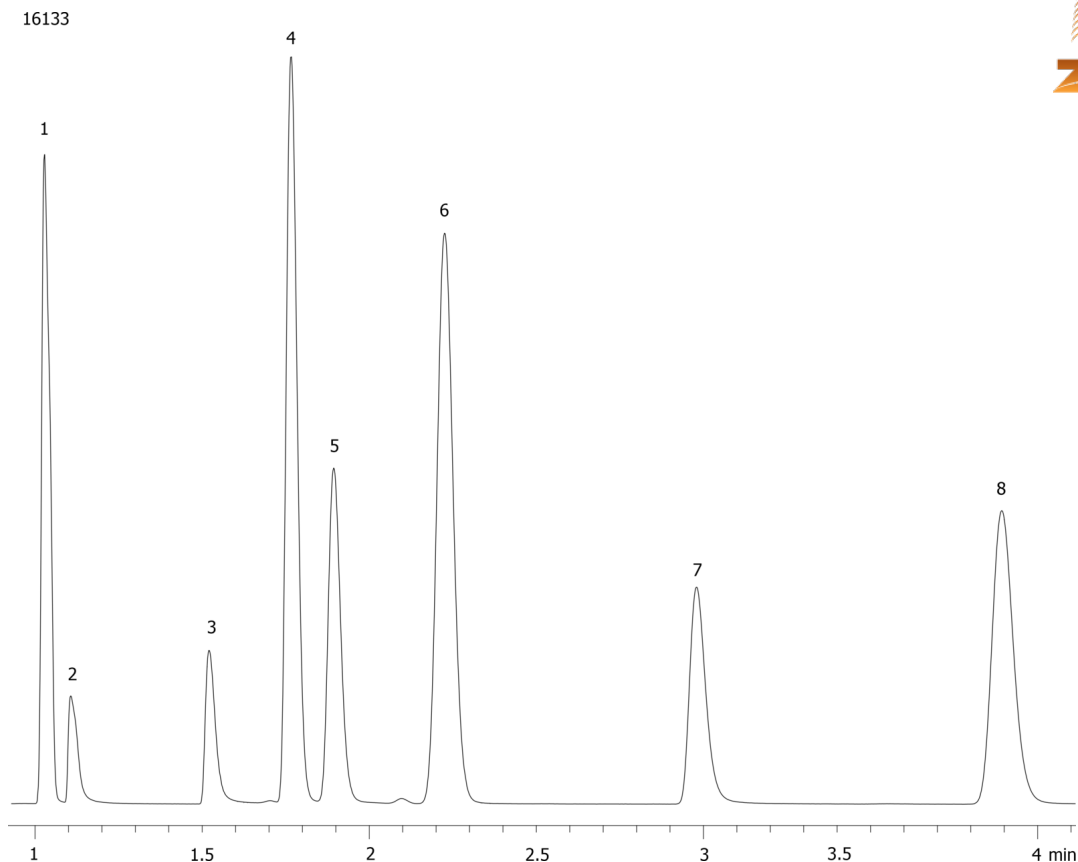