


Blood Alcohols on ZB-1701**Column:** ZebtronTM ZB-1701, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:** 14% Cyanopropylphenyl 86% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G006-11**Oven Profile:** 30°C Isothermal**Carrier Gas:** Constant Flow Helium, 1.8 mL/min**Injection:** Split 51:1 200 μ L @ 200°C**Detection:** Flame Ionization (FID) (220°C)**Analyst Note:** Static Headspace: 1mL 50ppm sample in 20mL vial at 70°C for 15min.**Products used in this application:****ANALYTES:**

- 1 Acetaldehyde
- 2 Methanol
- 3 Ethanol
- 4 Acetone
- 5 Isopropanol
- 6 t-Butanol
- 7 n-Propanol
- 8 2-Butanol

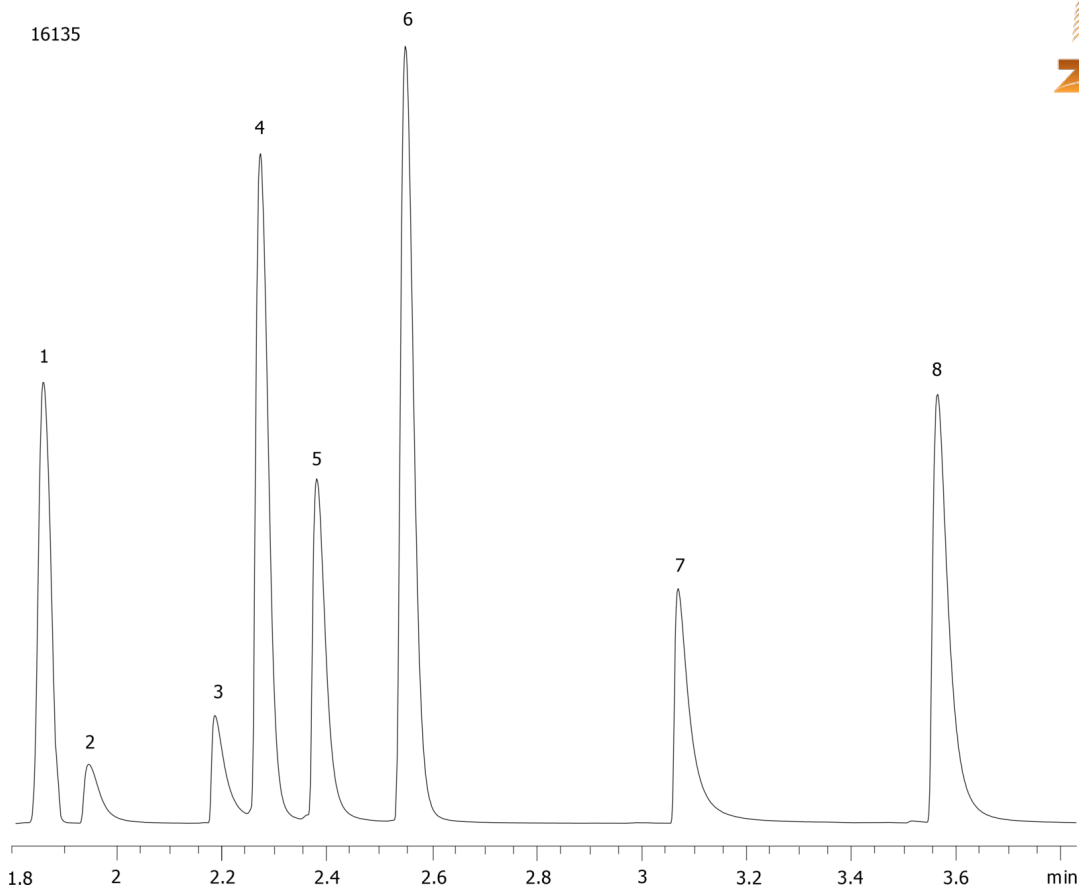