


DDD, DDT, And DDE Isomers on ZB-MR2**Column:** Zebron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:****Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:** 100 °C for 0.5 min to 200 °C at 25 °C/min to 320 °C at 15 °C/min for 2 min**Carrier Gas:** Constant Pressure Helium, 0.8 mL/min**Injection:** Splitless :1 1 μ L @ 260°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 25 ng on column**Products used in this application:****ANALYTES:**

- 1 o,p'-DDE
- 2 p,p'-DDE
- 3 o,p'-DDD
- 4 o,p'-DDT
- 5 p,p'-DDD
- 6 p,p'-DDT

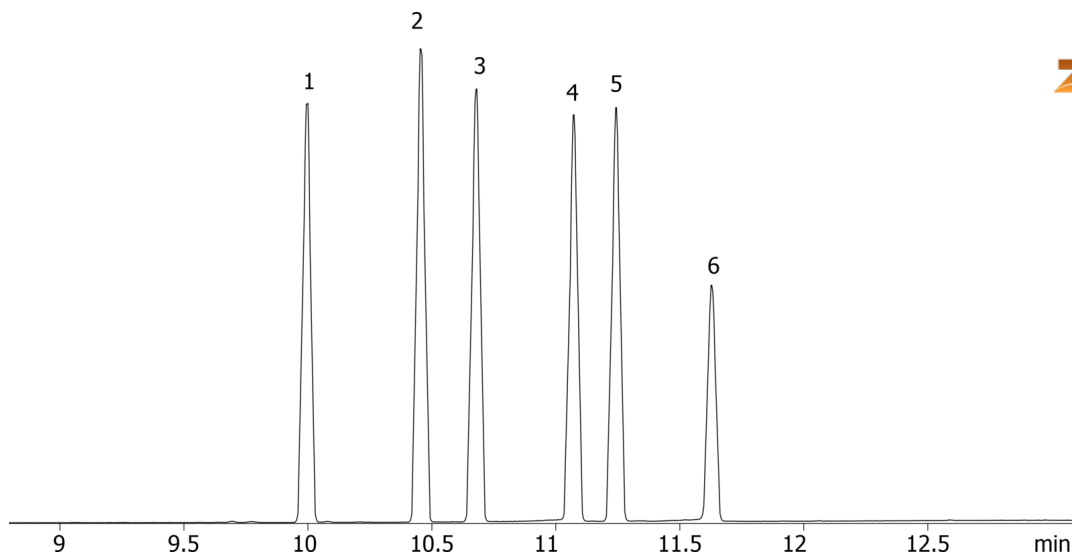