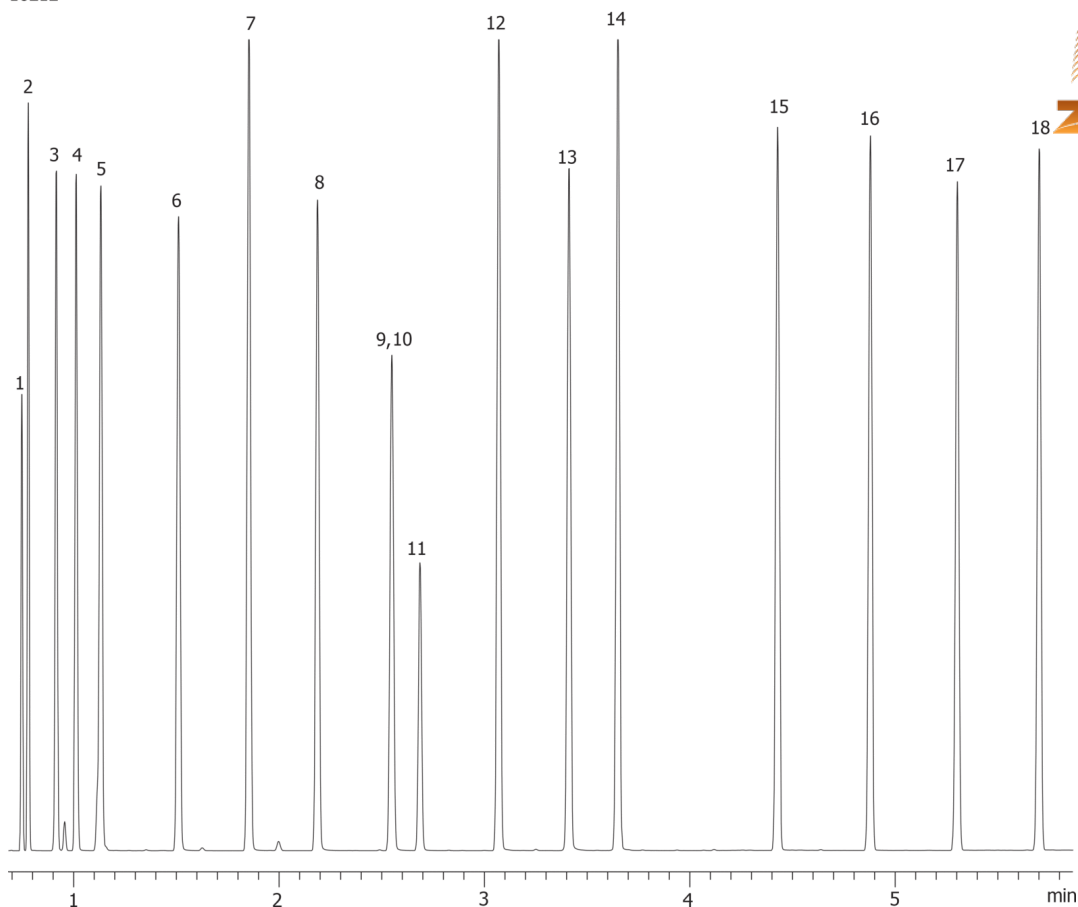


Fast ASTM D3710 on ZB-1**Column:** Zebron[™] ZB-1, GC Cap. Column 15 m x 0.53 mm x 0.50 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.53 mm x 0.5 μ m**Order No:** 7EK-G001-17**Oven Profile:** 35 °C for 1 min. to 200 °C at 35 °C/min. for 1 min.**Carrier Gas:** Constant Pressure Helium, 7 mL/min**Injection:** Split 50:1 0.1 μ L @ 200°C**Detection:** Flame Ionization (FID) (220°C)

16212

**Products used in this application:**

Fast ASTM D3710 on ZB-1

ANALYTES:

- 1 Isopentane
- 2 Pentane
- 3 2-Methylpentane
- 4 Hexane
- 5 2,4-Dimethylpentane
- 6 Heptane
- 7 Toluene
- 8 Octane
- 9 m-Xylene
- 10 p-Xylene
- 11 o-Xylene
- 12 Propylbenzene
- 13 Decane
- 14 Butylbenzene
- 15 Dodecane
- 16 Tridecane
- 17 Tetradecane
- 18 Pentadecane

