


Food and Flavors: Bergamot Oil on ZB-MR1**Column:** Zebtron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 40 °C for 3 min. to 240 °C at 5 °C/min. for 3 min.**Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 200:1 1 μ L @ 250°C**Detection:** Mass Spectrometer (MS) (300°C)**Analyst Note:** Bergamot oil 1 %**Products used in this application:**

Food and Flavors: Bergamot Oil on ZB-MR1

ANALYTES:

- 1 alpha-Thujene
- 2 alpha-Pinene
- 3 Sabinene
- 4 beta-Pinene
- 5 beta-Myrcene
- 6 Limonene
- 7 para-Cymene
- 8 gamma-Terpinene
- 9 Linalool
- 10 Linalyl acetate
- 11 Neryl acetate
- 12 Geranyl acetate

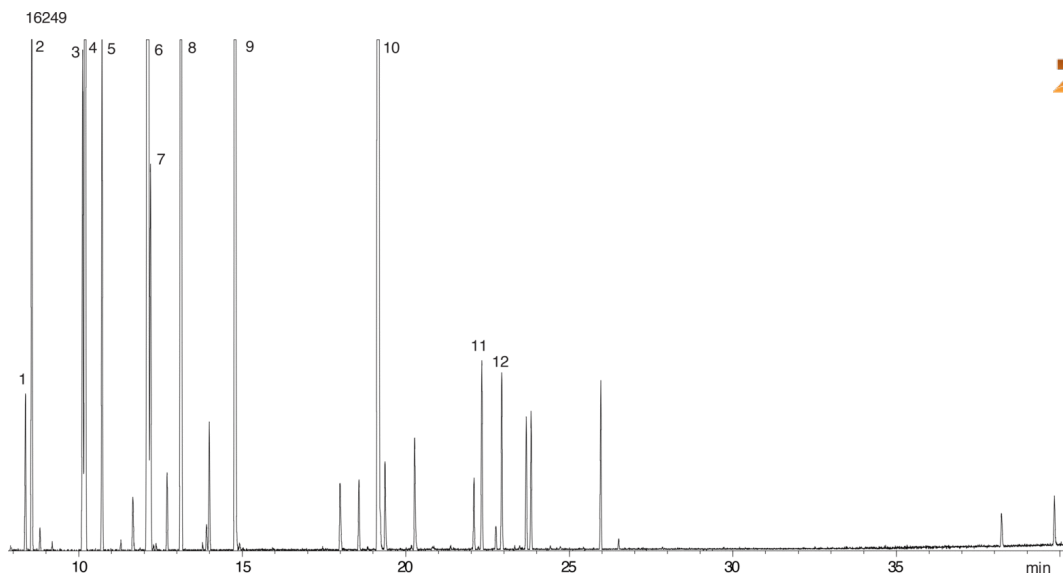