


ASTM Method D2887-97: Boiling Range Distribution of Petroleum Fractions on ZB-MR2**Column:** Zebron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:** 70 °C to 320 °C at 8 °C/min.**Carrier Gas:** Constant Flow Helium, 4 mL/min**Injection:** Split 10:1 0.5 μ L @ 310°C**Detection:** Flame Ionization (FID) (320°C)**Analyst Note:** Diesel Fuel 600ppm in Hexane**Products used in this application:**

ASTM Method D2887-97: Boiling Range Distribution of Petroleum Fractions on ZB-MR2

ANALYTES:

- 1 n-Nonane (C9)
- 2 n-Decane (C10)
- 3 n-Dodecane (C12)
- 4 n-Tetradecane (C14)
- 5 n-Hexadecane (C16)
- 6 n-Octadecane (C18)
- 7 n-Nonadecane (C19)
- 8 n-Eicosane (C20)
- 9 n-Docosane (C22)
- 10 n-Tetracosane (C24)
- 11 n-Hexacosane (C26)

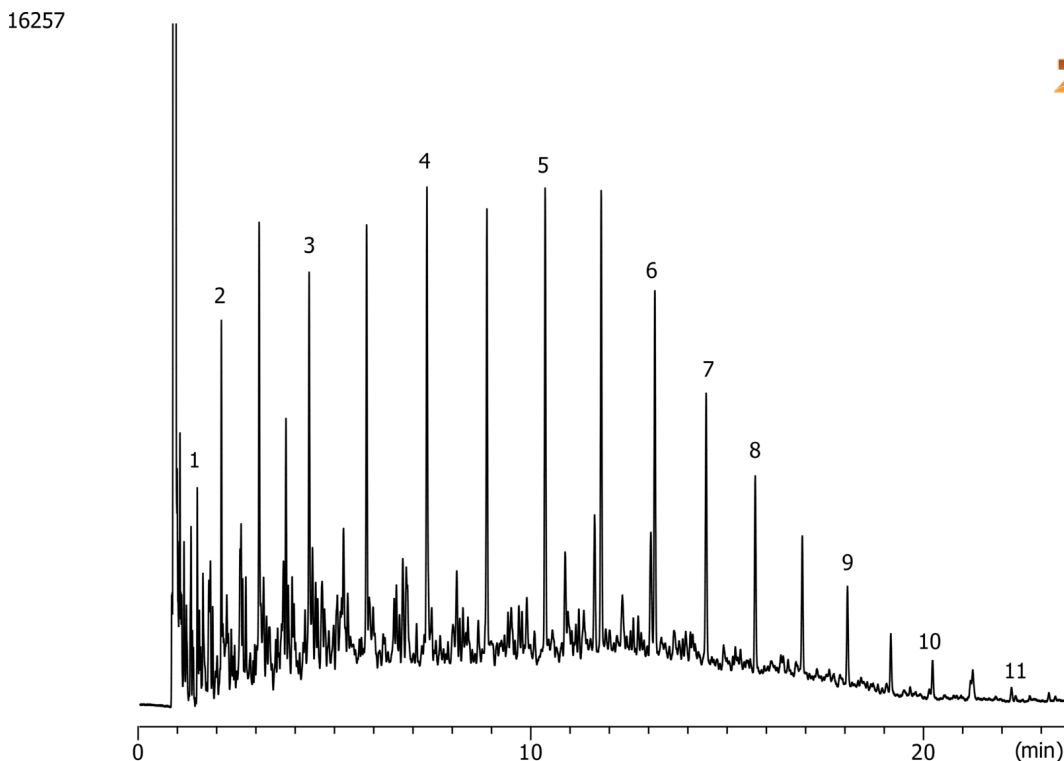