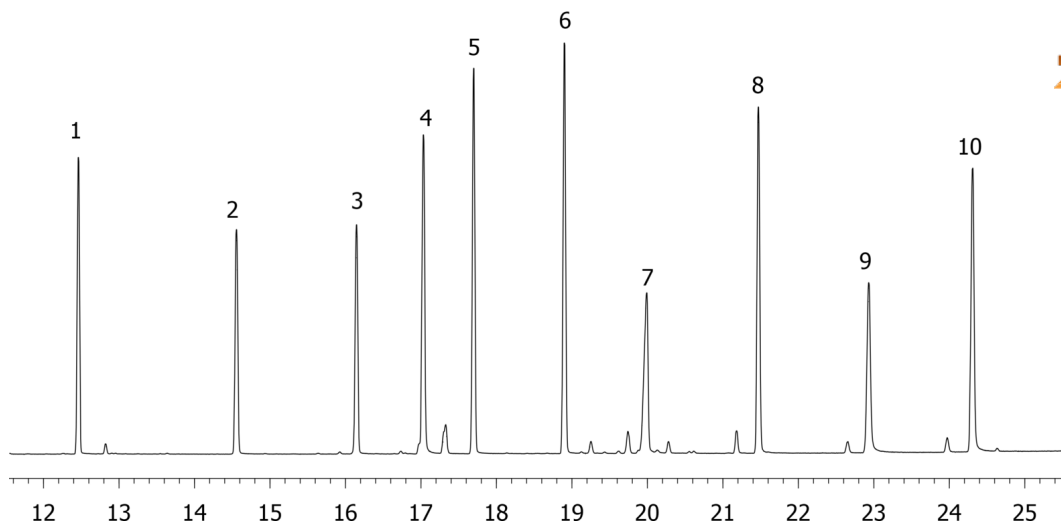


**Nitrogen and Phosphorous Pesticides on ZB-MR1****Column:** Zebtron<sup>™</sup> ZB-MultiResidue<sup>™</sup>-1, GC Cap. Column 30 m x 0.25 mm x 0.25  $\mu$ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.25  $\mu$ m**Order No:** 7HG-G016-11**Oven Profile:** 60 °C for 0.5 min to 320 °C @ 10 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Splitless :1 1  $\mu$ L @ 260°C**Detection:** Mass Selective (MSD) (180°C)

16375

**Products used in this application:**

## Nitrogen and Phosphorous Pesticides on ZB-MR1

### ANALYTES:

- 1 Butylate
- 2 Molinate
- 3 Chlorpropham
- 4 Atraton
- 5 Pronamide
- 6 Alachlor
- 7 Bromacil
- 8 Tetrachlorvinphos
- 9 Tricyclazole
- 10 Hexazinone

