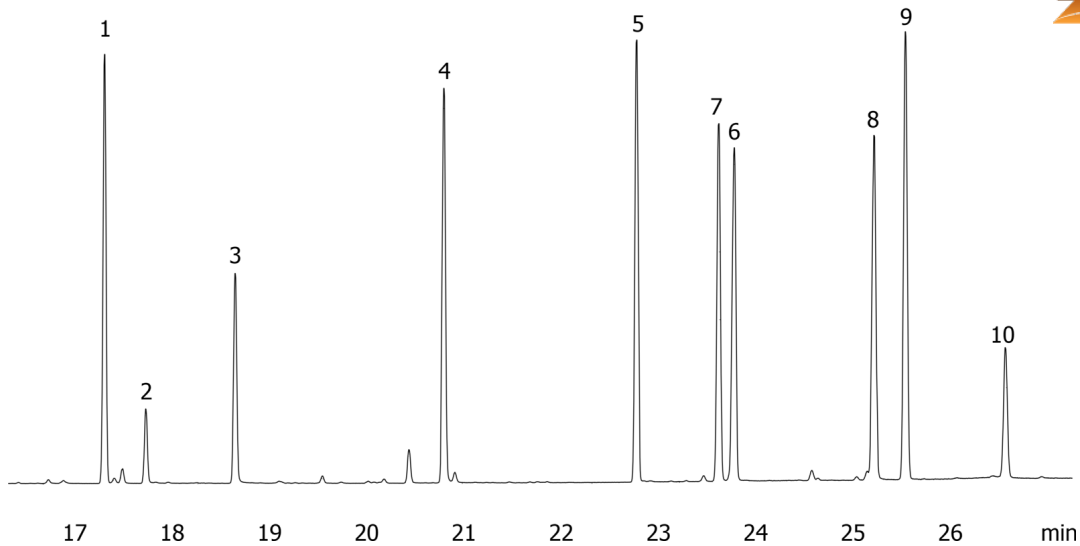


Nitrogen and Phosphorous Pesticides on ZB-MR1**Column:** Zebtron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 60 °C for 0.5 min to 320 °C @ 10 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Splitless :1 1 μ L @ 260°C**Detection:** Mass Selective (MSD) (180°C)

16378

**Products used in this application:**

Nitrogen and Phosphorous Pesticides on ZB-MR1

ANALYTES:

- 1 Terbufos
- 2 Phosphamidon
- 3 Phosphamidon
- 4 Chlorfenvinphos
- 5 Ethion
- 6 Famphur
- 7 Carbophenothion
- 8 Phosmet
- 9 Leptophos
- 10 Azinphos-ethyl

