


Nitrogen and Phosphorous Pesticides on ZB-MR2

Column: Zebtron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 µm, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.2 µm**Order No:** 7HG-G017-10**Oven Profile:** 60 °C for 0.5 min to 120 °C @ 10 °C/min to 240 °C @ 4 °C/min to 320 °C @ 15 °C/min for 1 min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Splitless :1 1 µL @ 260°C**Detection:** Mass Selective (MSD) (180°C)

16381



Products used in this application:

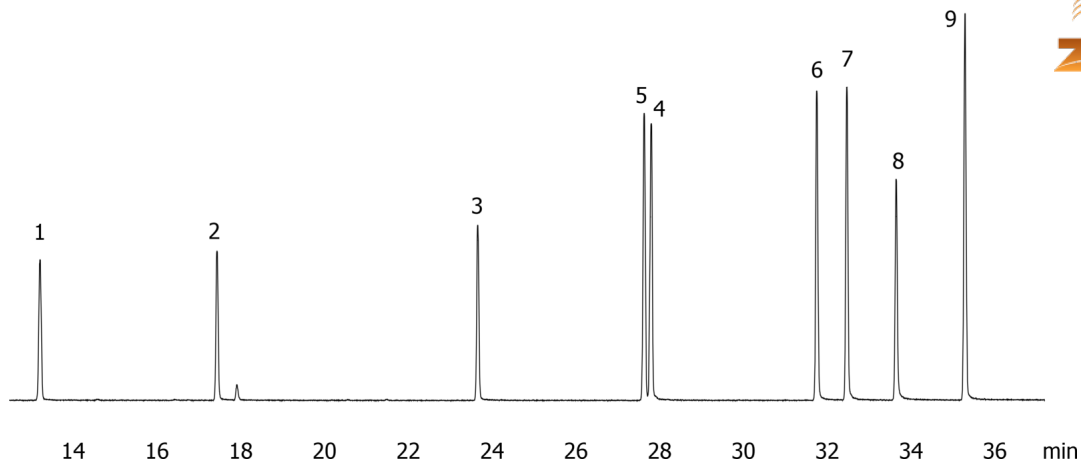

ANALYTES:

- 1 EPTC
- 2 Mevinphos
- 3 Ethoprop
- 4 Atrazine
- 5 Propazine
- 6 Prometryn
- 7 Terbutryn
- 8 Triadimefon
- 9 Diphenamide

