


Nitrogen and Phosphorous Pesticides on ZB-MR2**Column:** Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:** 60 °C for 0.5 min to 120 °C @ 10 °C/min to 240 °C @ 4 °C/min to 320 °C @ 15 °C/min for 1 min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Splitless :1 1 μ L @ 260°C**Detection:** Mass Selective (MSD) (180°C)

16384

**Products used in this application:****ANALYTES:**

- 1 DEET
- 2 Terbacil
- 3 Metribuzin
- 4 Bromacil
- 5 Triadimefon
- 6 Tricyclazole
- 7 Hexazinone

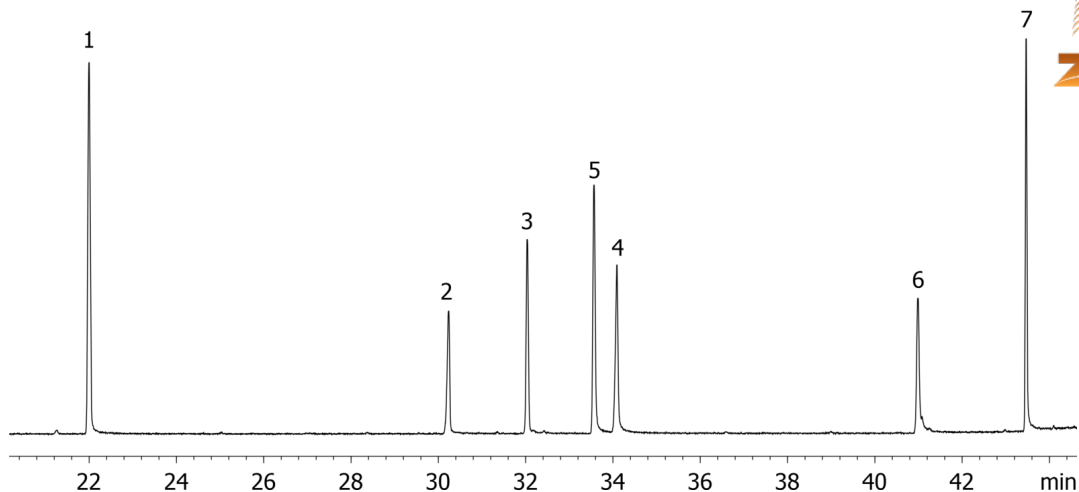