

Cytochrome C Mixture on Jupiter 3u C18 40°C

Column: Jupiter® 3 µm C18 300 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4263-E0

Elution Type: Gradient

Eluent A: 0.1% TFA and 5% Acetonitrile in Water

Eluent B: 0.1% TFA in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	15	45	55

Flow Rate: 1 mL/min


Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Diode Array (PDA) @ 214 nm (40 °C)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Widespore C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4321


- Holder Part No.: KJ0-4282

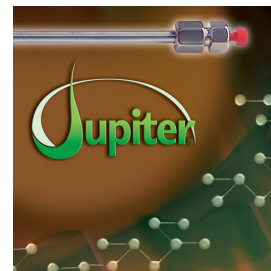

Products used in this application:



16451

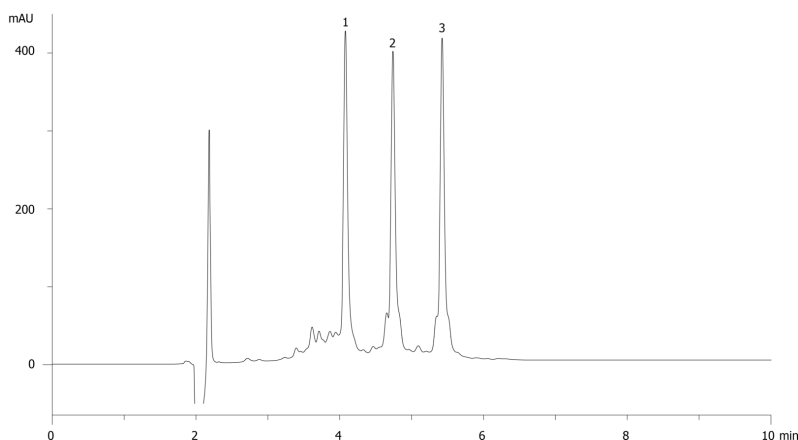

ANALYTES:

- 1** Cytochrome c Horse Heart
Retention Time: 4.08 min
- 2** Cytochrome c Bovine Heart
Retention Time: 4.74 min
- 3** Cytochrome c Saccharomyces Cerevisiae
Retention Time: 5.43 min

