


Chlorophenols on ZB-MR1

Column: Zebtron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G016-11**Oven Profile:** 60 °C for 0.5 min to 130 °C @ 10 °C/min to 210 °C @ 4 °C/min to 300 °C @ 10 °C/min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Splitless :1 1 µL @ 260°C**Detection:** Mass Selective (MSD) (180°C)**Analyst Note:** Column with 5 meter guard. Analytes are 0.5 ppm.
Application courtesy of Annegret Tischendorf, Südsachsen Wasser GmbH

Products used in this application:

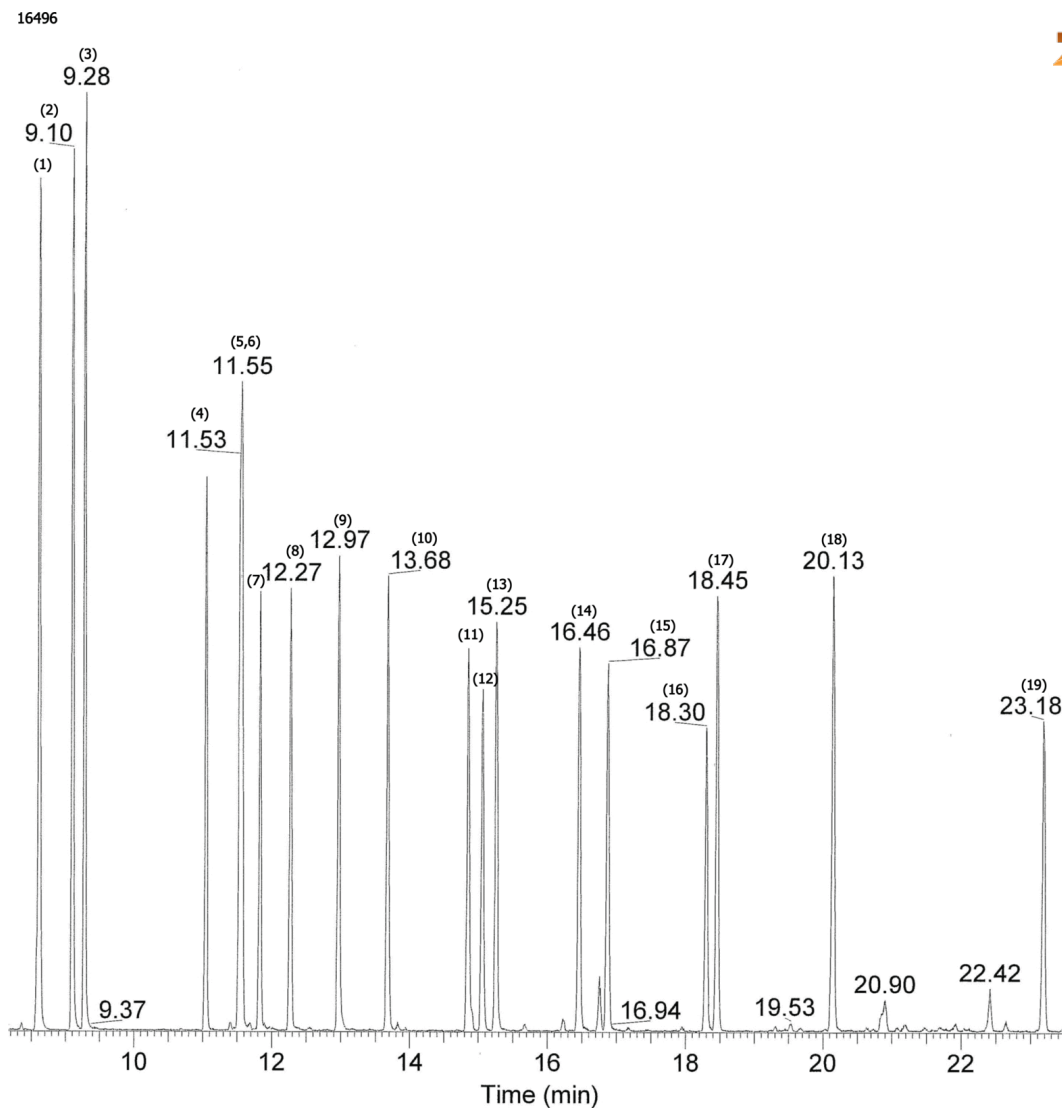

Chlorophenols on ZB-MR1

ANALYTES:

- 1** 2-Chlorophenol
- 2** 3-Chlorophenol
- 3** 4-Chlorophenol
- 4** 2,6-Dichlorophenol
- 5** 2,4-Dichlorophenol
- 6** 2,5-Dichlorophenol
- 7** 3,5-Dichlorophenol
- 8** 2,3-Dichlorophenol
- 9** 3,4-Dichlorophenol
- 10** 2,4,6-Trichlorophenol
- 11** 2,3,6-Trichlorophenol
- 12** 2,3,5-Trichlorophenol
- 13** 2,4,5-Trichlorophenol
- 14** 2,3,4-Trichlorophenol
- 15** 3,4,5-Trichlorophenol
- 16** 2,3,5,6-Tetrachlorophenol
- 17** 2,3,4,6-Tetrachlorophenol
- 18** 2,3,4,5-Tetrachlorophenol
- 19** Pentachlorophenol

