


Polynuclear Aromatic Hydrocarbons (PAHs) on ZB-MR1**Column:** Zebron™ ZB-MultiResidue™-1, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G016-11**Oven Profile:** 60 °C for 2 min to 200 °C @ 15 °C/min to 320 °C @ 30 °C/min for 2 min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Splitless :1 1 µL @ 260°C**Detection:** Mass Selective (MSD) (180°C)**Analyst Note:** Column with 5 meter guard. Analytes are 1.0 ppm.
Application courtesy of Annegret Tischendorf, Südsachsen Wasser GmbH**Products used in this application:**

Polynuclear Aromatic Hydrocarbons (PAHs) on ZB-MR1

ANALYTES:

- 1 Naphthalene
- 2 Acenaphthalene
- 3 Acenaphthene
- 4 Fluorene
- 5 Phenanthrene
- 6 Anthracene
- 7 Fluoranthene
- 8 Pyrene
- 9 Benz[a]anthracene
- 10 Chrysene
- 11 Benzo[b]fluroanthene
- 12 Benzo[k]fluroanthene
- 13 Benzo[a]pyrene
- 14 Dibenzo[a,h]anthracene
- 15 Indeno[1,2,3-cd]pyrene
- 16 Benzo[g,h,i]perylene

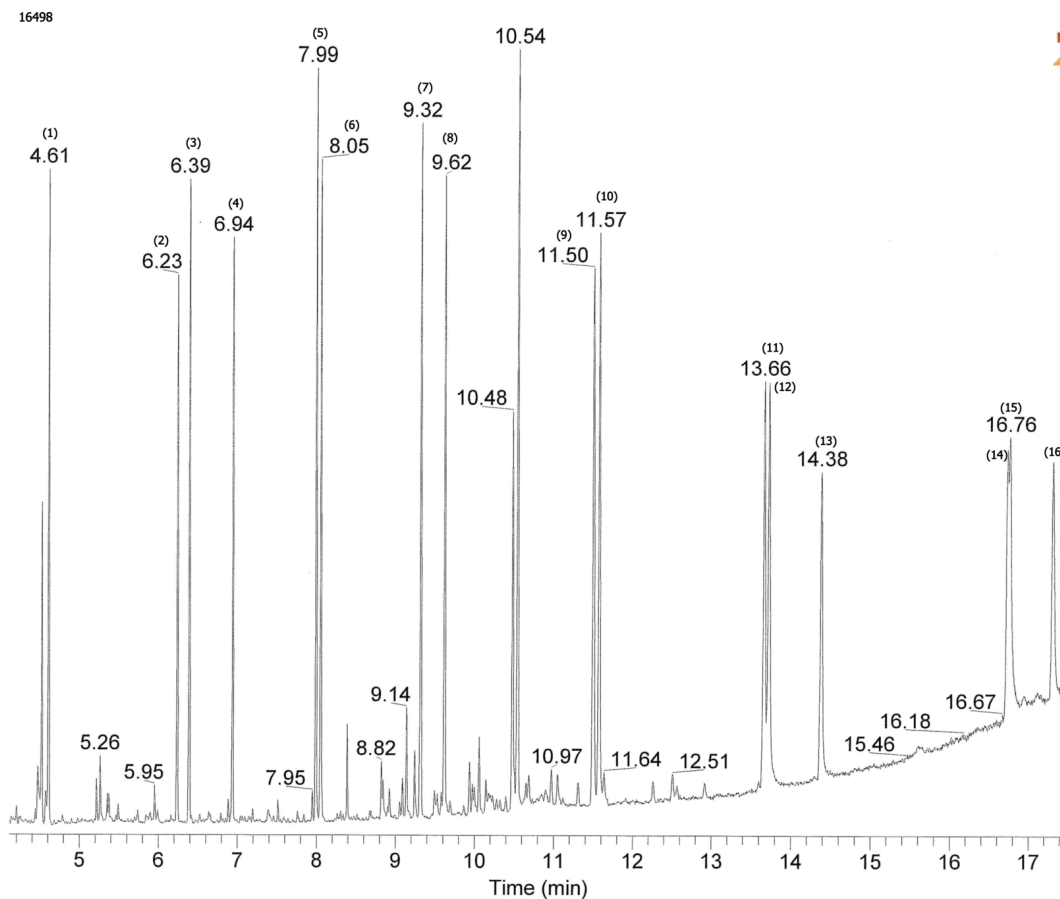